

# Joint British Society consensus recommendations for magnetic resonance imaging for patients with cardiac implantable electronic devices

Anish Bhuvu <sup>1,2</sup>, Geoff Charles-Edwards,<sup>3,4</sup> Jonathan Ashmore,<sup>5,6</sup> Alexandra Lipton,<sup>7</sup> Matthew Benbow,<sup>8,9</sup> David Grainger,<sup>10</sup> Trudie Lobban,<sup>11</sup> Deepa Gopalan,<sup>12,13</sup> Alistair Slade,<sup>14,15</sup> Giles Roditi,<sup>16,17</sup> Charlotte Manisty <sup>18,19,20</sup>

For numbered affiliations see end of article.

## Correspondence to

Dr Charlotte Manisty, Institute of Cardiovascular Science, University College London, London, UK; c.manisty@ucl.ac.uk

## ABSTRACT

Magnetic Resonance Imaging (MRI) is increasingly a fundamental component of the diagnostic pathway across a range of conditions. Historically, the presence of a cardiac implantable electronic device (CIED) has been a contraindication for MRI, however, development of *MR Conditional* devices that can be scanned under strict protocols has facilitated the provision of MRI for patients. Additionally, there is growing safety data to support MR scanning in patients with CIEDs that do not have MR safety labelling or with *MR Conditional* CIEDs where certain conditions are not met, where the clinical justification is robust. This means that almost all patients with cardiac devices should now have the same access to MRI scanning in the National Health Service as the general population. Provision of MRI to patients with CIED, however, remains limited in the UK, with only half of units accepting scan requests even for patients with *MR Conditional* CIEDs. Service delivery requires specialist equipment and robust protocols to ensure patient safety and facilitate workflows, meanwhile demanding collaboration between healthcare professionals across many disciplines. This document provides consensus recommendations from across the relevant stakeholder professional bodies and patient groups to encourage provision of safe MRI for patients with CIEDs.

## SCOPE

The aim of this joint multiprofessional societal guidance is to provide consensus recommendations for best practice management of patients with cardiac implantable electronic devices (CIEDs), who require investigation using MRI in the UK. With representation from all involved in the patient pathway (including patients), we aim to highlight areas of clear recommendations which should be adhered to, alongside consensus recommendations where no current guidelines exist or are perceived to conflict. This document provides a recommended protocol and workflow alongside specific guidance for the different personnel involved in the clinical pathway, outlining relevant responsibilities and procedures. Additionally, the risks associated with scenarios where particular conditions of *MR Conditional* CIEDs are not met or where the CIED system does not currently have regulatory approval to undergo MRI are summarised, to aid a local decision to scan patients in these scenarios where

the clinical benefit outweighs the risk (figure 1). Suggested statements for consent of these patients are also provided. For each section throughout this document, points are presented as statements with distinction made between mandated and consensus recommendations.

This guidance aims to support the development of new providers of MRI services to patients with CIED (adults and children) and help the growth of existing services, in order to facilitate equitable provision for those patients in need wherever they may be in the UK. Guidance on MRI scanning of patients with CIEDs from other professional bodies have been considered in the production of these UK recommendations.<sup>1–8</sup> This guidance is not intended to provide a comprehensive literature review, which can be found elsewhere.<sup>1,9</sup>

## BACKGROUND


Between 2016 and 2017, 3.7 million MRI scans were performed in England.<sup>10</sup> MRI is one of the fastest growing imaging modalities with many diagnostic and treatment pathways increasingly dependent on MRI, including orthopaedics, neurosurgery and radiotherapy.<sup>10</sup> Alongside this, implantation rates of CIEDs are rising—there are currently half a million people in the UK with cardiac permanent pacemakers or implantable cardioverter-defibrillators (ICDs) and over 40 000 new implants per year.<sup>11</sup> These patients have historically been prevented from having MRI scans because of safety concerns, although half of this patient group are aged over 65 years and therefore have high clinical requirement for imaging due to comorbidities.<sup>12</sup> Data suggest that 7%–17% of patients undergoing device implantation have MRI scans requested in the first 12 months post device implantation, highlighting the imperative to enable scanning where feasible.<sup>13,14</sup> This demand for MRI in patients with CIEDs is growing rapidly at an estimated 10 000 scans a year based on annual growth in CIED implantations and MRI requests,<sup>10,11,15</sup> highlighting the requirement for consensus recommendations for service delivery.

In response to this, industry has adapted the hardware and software in CIEDs to develop *MR Conditional* devices with regulatory approval for MRI scanning under strict conditions. Almost all CIED



© Author(s) (or their employer(s)) 2022. No commercial re-use. See rights and permissions. Published by BMJ.

**To cite:** Bhuvu A, Charles-Edwards G, Ashmore J, et al. *Heart* Epub ahead of print: [please include Day Month Year]. doi:10.1136/heartjnl-2022-320810

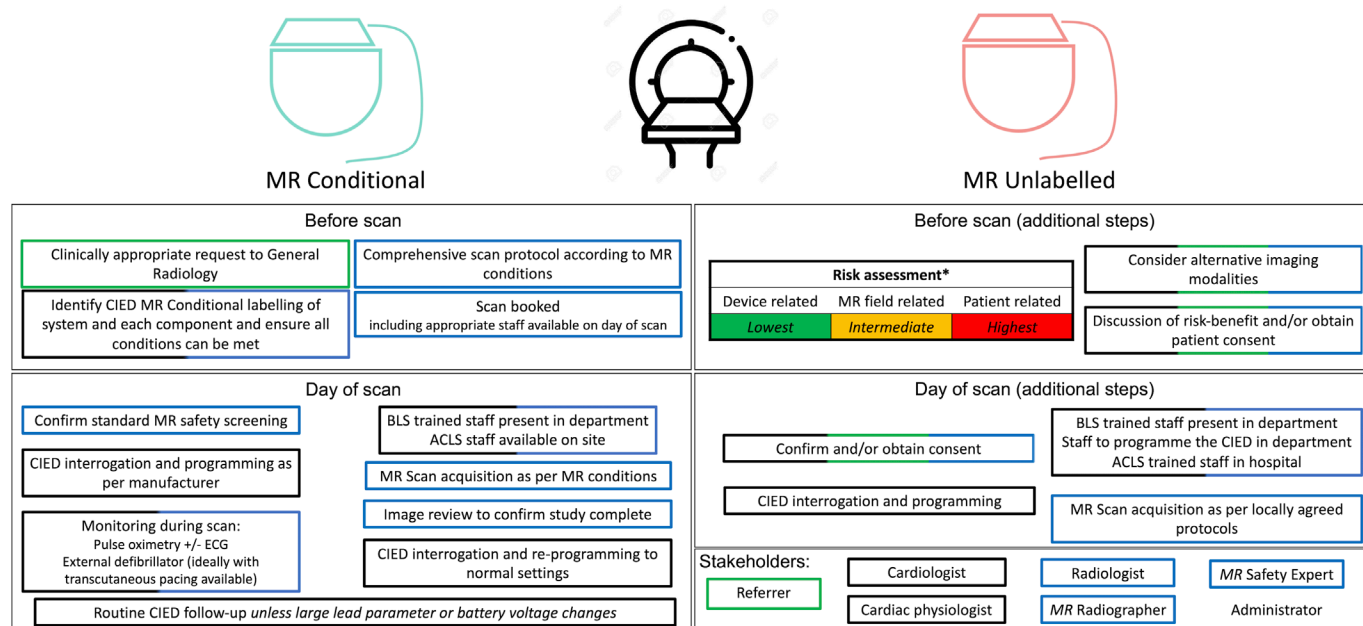

**Figure 1** Recommendations for MRI workflows in patients with cardiac implantable electronic devices. \*Higher risk scenarios include the presence of fractured, epicardial, abandoned leads; recent implantation; battery at elective replacement indicator; deactivated systems; lead parameters outside manufacturer recommendations, other implants present. ACLS, adult cardiac life support; BLS, basic life support; SAR, specific absorption rate.

types are now available in *MR Conditional* models, and manufacturers report that currently almost all new CIED implants in the UK are *MR Conditional*. Provided all MR conditions are met including device reprogramming, these patients can safely undergo MR scanning. Alongside this, there is increasing clinical evidence that the risks associated with scanning a patient with a CIED that has not been formally tested and approved to undergo MRI, or where certain conditions of *MR Conditional* CIEDs are unmet, are lower than previously thought provided scans are performed under similar strict conditions to those required for *MR Conditional* CIEDs. Importantly, where there are no appropriate alternatives, MRI scanning is commonly justified and should be considered.

Barriers exist however at multiple levels from referrer to reporting radiologist—patients with CIEDs are approximately 50 times less likely to be referred for MRI than the general population, and workflows need to incorporate time and collaboration from multiple hospital departments with no established funding strategies that recognise service complexity.<sup>15–18</sup> Progress has been made—a joint statement by the Clinical Imaging Board and British Cardiovascular Society demonstrates high-level consensus that new working practices are required.<sup>19</sup> The 2018 British Heart Rhythm Society Standards for Implantation and Follow-up of Cardiac Rhythm Devices explicitly require CIED implantation centres to provide pacing support for MRI units.<sup>20</sup> A recent survey of MRI departments in England showed that challenges to provision of MRI to patients with CIEDs persist—only 53% of units will scan patients with *MR Conditional* devices, and there remains an estimated 10-fold service underprovision.<sup>15 21</sup>

### Clinical demand for MRI: a changing landscape

MRI has evolved as a powerful and versatile diagnostic imaging modality since its introduction into clinical use in the early 1980s. Technological advances have led to clinical application for diagnosis and treatment planning across all body areas and systems with profound impact on patient care. Consistent growth in

referrals for MRI reflects the expanding clinical indications and incorporation into many guideline-recommended clinical diagnostic pathways.<sup>10</sup> In the acute setting, timely provision of MR imaging is fundamental for diagnosis of a variety of conditions including acute ischaemic brain injury, spinal cord compression, spinal infection and trauma, while MRI is increasingly indicated in the oncology setting where it is first line for detection, characterisation and staging of many tumours including suspected clinically localised prostate cancer.<sup>22</sup>

Where MRI is unavailable, clinicians are generally forced to opt for alternate investigations that may be more invasive or have lower diagnostic accuracy, resulting in late or misdiagnosis with inherent clinical complications. Similarly, treatments that require pretherapy MRI planning such as neurosurgical procedures or MRI-guided stereotactic radiotherapy are unavailable to those patients where MRI is contraindicated, potentially impacting clinical outcomes.

### TERMINOLOGY

The following terms are defined by the international standard ASTM F2503-20, and are recognised and used globally by the MR community, medical device manufacturers and regulatory bodies.

- **MR Environment**—the three-dimensional volume of space surrounding the MR magnet that contains both the Faraday shielded volume and the 0.50 mT field contour (5 G line). This volume is the area in which an item might pose a hazard from exposure to the electromagnetic fields produced by the MR equipment and accessories.
- **MR Safe**—an item that poses no known hazards resulting from exposure to any *MR Environment*. *MR Safe* items are composed of materials that are electrically non-conductive, non-metallic and non-magnetic.
- **MR Conditional**—an item with demonstrated safety in the *MR Environment* within defined conditions including conditions for the static magnetic field, the time-varying gradient magnetic fields and the radiofrequency (RF) fields.

No CIEDs are *MR Safe*, since they contain materials that are electrically conductive. However, the vast majority of CIEDs now implanted in the UK are now *MR Conditional*. This means that if all of the stated MR conditions are met, the manufacturer of the CIED is providing assurance that in terms of the MR safety issues related to that specific device, it is safe for the patient to undergo MRI.

- ▶ *MR Unsafe*—an item which poses unacceptable risks to the patient, medical staff or other persons within the *MR Environment*.

The Medicines and Healthcare Products Regulatory Agency (MHRA) guidelines for MR safety,<sup>23</sup> the primary reference for MR safety guidance in the UK, additionally define the following term:

- ▶ *MR Unlabelled*—an item without an *MR Safe*, *MR Conditional* or *MR Unsafe* label.

In the context of CIEDs, *MR Unlabelled* items have also been described as 'conventional', 'legacy', 'MR non-conditional' and 'non-MR Conditional'.

Some *MR Unlabelled* items will clearly be unsafe in the *MR Environment*, for example, a ferromagnetic gas cylinder. Others will clearly be safe, for example, a saline bag. For many *MR Unlabelled* items, the MR safety risks will lie somewhere between these two extremes and may not be fully understood, particularly for implants. A key aspect of *MR Unlabelled* items is that no statement about MR safety is being made by the manufacturer of the item. Consequently, a local decision is required on whether to bring such items into the *MR Environment* based on a risk-benefit assessment.

Importantly, for *MR Conditional* items in scenarios where any of the MR conditions are unmet, again no statement about MR safety is being made. In such scenarios, these items should be managed in the same way as *MR Unlabelled* items, requiring a local decision to be made based on a risk-benefit assessment.<sup>23</sup> Regarding CIEDs, the device system needs to be assessed in its entirety—as a whole—generator, lead(s) and any other system component in combination. Although individual components may be *MR Conditional*, manufacturers clearly stipulate the generator-lead combinations that have been tested and approved to be *MR Conditional*, and devices outside of these recommendations should be considered *MR Unlabelled*.

The MHRA safety guidelines defines the following terms that are used in this guidance<sup>23</sup>:

- ▶ *MR Responsible Person*—someone who takes on the day-to-day responsibility for MR safety.
- ▶ *MR Safety Expert*—someone who can adequately advise on the necessary engineering, scientific and administrative aspects of MR safety. Their knowledge of MR physics should enable them to advise on the risks associated with individual procedures and on methods to mitigate these risks.

For the purpose of this guidance, some of the tasks may be undertaken by persons other than the MR Safety Expert but who have the required scientific knowledge.

- ▶ *MR Operator*—someone who is entitled to operate the MRI equipment. MR Operators are typically MR Radiographers, but may be assistant practitioners, radiologists, cardiologists or physicists.

For the purpose of this guidance, the term MR Radiographer will be used as it is recognised in the UK that Health and Care Professions Council (HCPC) registered radiographers perform the overwhelming majority of diagnostic MR scans. Where the MR Operator is not an HCPC registered radiographer services must ensure that clear governance processes are in place

outlining roles, scope, supervision and responsibilities including responsibility for the safety of the patient during scanning.

Finally, the following term is defined in this guidance:

- ▶ *MR Clinician*—any clinician responsible for reviewing appropriateness of referrals, protocolling and/or reporting the MRI scan of a patient with a CIED. For most sites these will be radiologists, but this may vary dependent on scan indication and setting, for example, a cardiologist for a cardiology-led MRI service or appropriately trained reporting radiographers. Departments scanning cardiac devices should ensure that more than one MR Clinician is trained and familiar with processes and procedures required for the safety and workflows of cardiac device MRI.

## MR CONDITIONAL CIEDS


Historically, CIEDs were viewed as an absolute contraindication for MRI due to their perceived sensitivity to the strong static and time-varying magnetic fields produced by MRI scanners. These fields interact with medical devices in multiple ways, giving rise to various risks including mechanical forces (attraction, torque, vibration), heating, unintended stimulation and device malfunction. A number of technical developments have been incorporated into CIED design to mitigate these risks, including a reduction in the amount of ferromagnetic material, improved lead design and adapted software programming modes.<sup>24</sup> In Europe, this resulted in the approval of *MR Conditional* pacemakers in 2008 and *MR Conditional* ICDs in 2014, with subsequent introductions into North America a few years later. Modern implantable cardiac monitors (including implantable loop recorders and implantable pulmonary artery pressure monitors) are *MR Conditional*, and this has been the case for the commonly implanted models for over a decade. Importantly, the conditions associated with *MR Conditional* cardiac monitors are relatively simple to meet without the need for the cardiology support that is required for *MR Conditional* pacemakers and ICDs, and many devices currently do not require data download prior to scanning.

Provided all the MR conditions are met, *MR Conditional* devices have been demonstrated to be safe for patients to undergo MR scanning and have regulatory approval as such. Various studies have demonstrated no clinically significant complications in patients with *MR Conditional* CIEDs randomised to MRI.<sup>25–28</sup> Since their general introduction to our knowledge, there have been no adverse incidents associated with *MR Conditional* CIEDs undergoing MRI scanning when the MR conditions have been followed as per manufacturer recommendations.

## Recommendations: standardised protocol for all MR Conditional CIEDs

Manufacturer conditions for patients with *MR Conditional* CIEDs to undergo scans include considerations at the time of device implantation, scan booking, on the day of the scan prior to MRI, during the scan and after completing the study (figure 2). These recommendations should be considered in addition to any MR conditions stated by the CIED manufacturer (routinely available on manufacturer websites). Services should consider how best to setup their systems to manage time points such that patient are triaged safely and effectively. In particular, attention should be paid to who is authorising, aware and facilitating these patients being booked.



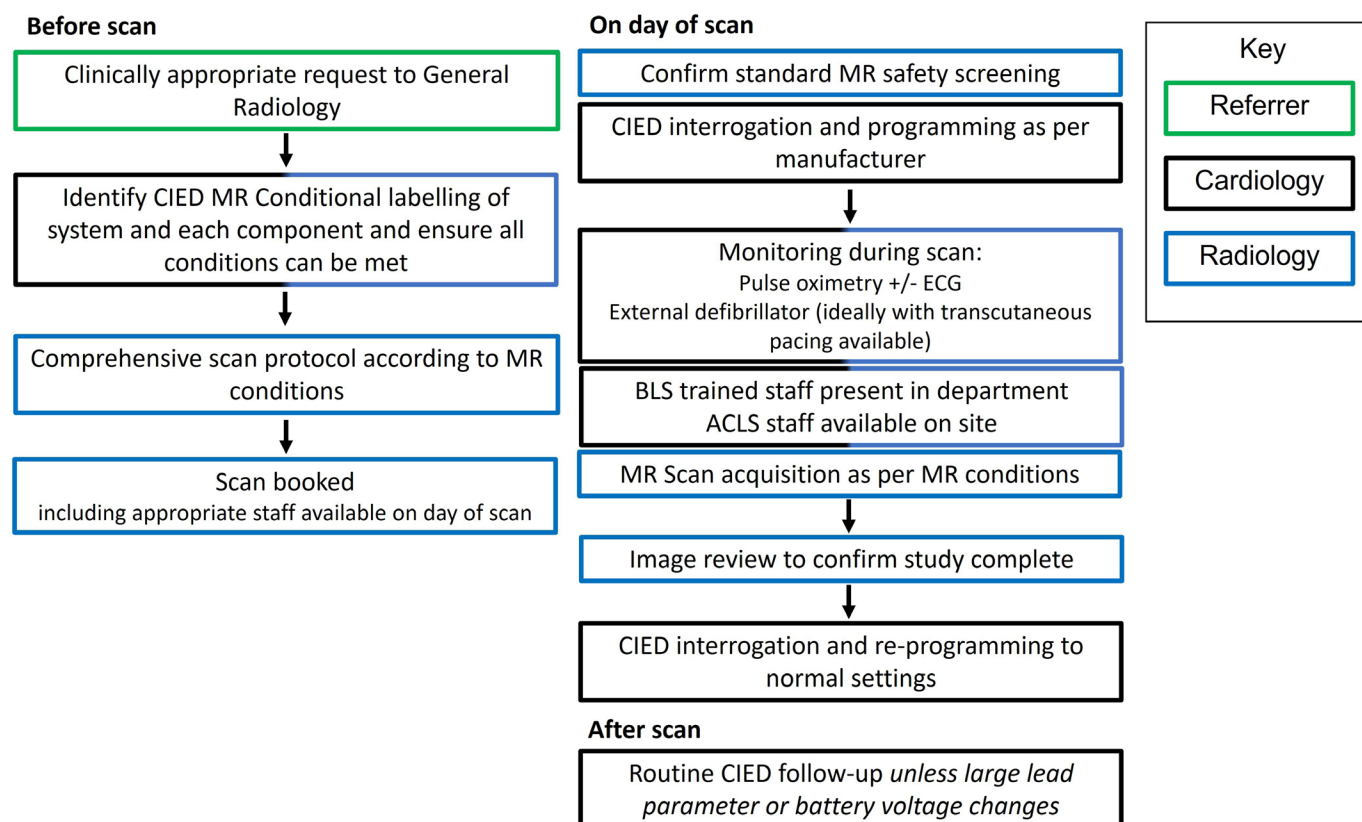

**Figure 2** Workflow for provision of MRI to patients with MR Conditional cardiac implantable electronic devices. ACLS, adult cardiac life support; BLS, basic life support; CIED, cardiac implantable electronic device; SAR, specific absorption rate.

#### Requirements at time of CIED implantation

All CIED implanting hospitals must adhere to British Heart Rhythm Society (BHRS) Standards for Device Implanting Centres to facilitate equitable provision of MRI to their patients.<sup>20</sup> It is recommended that *MR Conditional* CIED systems (generator and leads in combination) are the default selection for all new implantations, unless there is strong reason to do otherwise, or there are other absolute contraindications to MR scanning. Similarly, at the time of upgrade or generator change, implanting cardiologists should ensure that lead and generator combinations maintained to be from the same manufacturer to ensure that conditions of the *MR Conditional* CIED can still be met. The basic implications of an *MR Conditional* CIED should be explained to the patient at the time of implantation, and there should be written confirmation on the device identification card

stating whether the implanted CIED system is *MR Conditional* or not. This information should also be accessible in the medical notes of the patient. Each device has additional conditions that need to be fulfilled before a scan be performed (current device parameters acceptable, scanner and protocol conditions, etc), however these would be determined at the time of scan for patients with *MR Conditional* CIEDs.

For the rare scenarios where a fully *MR Conditional* CIED cannot be implanted (eg, where a patient's cardiac anatomy necessitates a specific lead choice leading to manufacturer mismatch between the leads and generator, or the device must be implanted outside of the pectoral region), discussion with the patient regarding the risks and benefits should take place prior to device implantation and included in the formal consent process by the implanting cardiologist. Suggested consent statements are

**Table 1** Suggested consent statement for implanting cardiologists when potentially implanting an *MR Unlabelled* device, or revising/upgrading a device that may then become MR Unlabelled (such as upgrading a pacemaker to a CRT using a lead from a different manufacturer to the existing leads and generator)

Implanting an <i>MR Unlabelled</i> CIED or revising an existing device such that the MR conditions cannot be met.	<p>There may be a need to implant or upgrade your pacemaker/defibrillator with a device system that has not been formally approved to undergo MRI scanning by the manufacturer.</p> <p>This means that it may be more difficult for you to have an MRI scan in the future should you need one. Although almost all devices can be scanned, these are generally only done in specialist centres.</p> <p>A decision to go ahead with a MRI scan may be made after discussing the possible benefits, risks and alternatives with your referring doctor at such a time an MRI is requested.</p> <p>Serious complications related to MRI occur in &lt;1 in every 2000 patients (about 0.05%) with these devices if there are no other high-risk features. These include, but are not limited to:</p> <ul style="list-style-type: none"> <li>▶ damage to the cardiac device;</li> <li>▶ abnormal heart rhythms;</li> <li>▶ excessive tissue heating.</li> </ul>
---	---

CIED, cardiac implantable electronic device; CRT, cardiac resynchronisation therapy.

included in [table 1](#). It should be explained to patients that having generator-lead manufacturer mismatch is no longer an absolute contraindication for MRI if there is sufficient clinical justification for the scan, although access to scans is likely to be more challenging as they are likely only be undertaken in centres that scan *MR Unlabelled* CIEDs. Transvenous CIED lead extraction has a procedure-related major complication rate reported of 0.19%–1.8%, including a mortality of 0.19%–1.2%.<sup>29</sup> The risk of MRI scanning when appropriate protocols for *MR Unlabelled* CIEDs are adhered to appears to be very significantly lower than the risk of lead extraction, and so we would not recommend lead extraction solely to facilitate MRI.<sup>30–32</sup>

Each device centre must ensure that they have arrangements in place that allow patient access to MRI scanning. As per BHRS standards, CIED component details and *MR Conditional* labelling should be provided by the cardiology or device clinic to MRI departments on request to support the scanning process, and should be made available on easily accessible electronic healthcare record systems. For patients with *MR Unlabelled* CIEDs, arrangements may be required for MRI referrals to an external centre. Details of a registered referrals network can be found elsewhere ([www.mrimypacemaker.com](http://www.mrimypacemaker.com)).<sup>33</sup>

## Referral requirements and workflows

### Identifying referrals

For all patients with CIEDs, the presence of a device must be flagged in the referral under all circumstances. All MRI units should have an established process in place to accept MRI referrals for patients with CIEDs, or have an established relationship with an external centre with a pathway for referral if there is no access to cardiology support locally. Clinically appropriate scan requests for patients with *MR Conditional* devices should not be declined without clear advice regarding how to access agreed referral pathways for alternative external sites, and MRI units within an institution that has CIED implantation services should not refer externally for MRI for *MR Conditional* devices. The referral process should be easily accessible to referring clinicians, although individual hospitals are likely to develop their own local protocols for accepting referrals. As many current electronic systems block requests if fields are checked for the presence of a CIED, a separate standardised booking proforma may be required. Work is currently in progress to update Royal College of Radiologists recommendations for electronic requesting systems (Order Communications) to facilitate MRI requests for patients with CIEDs, including fields for communicating the device details (including manufacturer and model) and *MR Conditional* labelling. An appropriate subspecialty MR radiologist should be available for discussion of the risk-benefit and potential for imaging using alternative modalities for referrals of patients with *MR Unlabelled* CIEDs. Staff availability should be checked prior to scan booking.

### Device information to determine MR Conditional status

Device information should be provided to MRI units prior to booking the scan. This should include manufacturer, model and implantation date for generator and each of the leads. We recommend that the responsibility for obtaining the information lies principally with the referrer. MRI departments should not be obligated to accept referrals until they are satisfied that the data provided are sufficient to allow them to identify that the CIED is *MR Conditional* and to check the specific MR conditions for that device.

Ideally, device information should be recorded in the patient's electronic patient record. However, many patients will have scans requested at hospital sites remote from where their device is followed up. In this situation, device information is best obtained from the patient's usual CIED clinic, who should provide a copy of the implant report and/or the last device check. If this is not available, patients may provide a copy of their device identification card (issued at the time of implant), but the cardiology team will also need to ensure that all conditions are met at the time of interrogation prescan. All patients with CIEDs should be screened for the existence of abandoned leads or implantation of other metallic devices that may not result in MR conditions for the device being met. Where there is uncertainty related to the presence of additional hardware, a chest radiograph can be performed as part of screening, although this should not be required routinely.

An *MR Conditional* CIED must have both the leads and generators implanted by the same manufacturer in a combination that has been tested and verified to be safe within an *MR Environment*. An *MR Conditional* generator and *MR Conditional* leads from different manufacturers does not constitute an *MR Conditional* system. While the risk profile of such a combination may be comparable to a *MR Conditional* CIED,<sup>34</sup> scanning should currently be performed as per protocols for *MR Unlabelled* CIEDs, see section *MR Unlabelled* CIED systems and *MR Conditional* CIEDs not fulfilling specified conditions. Device manufacturers offer verification tools which are centralised online ([www.mrimypacemaker.com](http://www.mrimypacemaker.com)),<sup>33</sup> and manufacturer representatives can also assist in the process. Verification of *MR Conditional* labelling should be performed by the most appropriate member of the team, but may be the cardiac physiologist, cardiologist, MR Radiographer, radiologist or *MR Safety Expert*. Ultimately, the MR Radiographer is responsible for the safety of the patient they are scanning and for ensuring all local procedures have been completed before the patient is brought into the *MR Environment*.

## Prescan preparation

### Consent

Patients with *MR Conditional* CIEDs do not require written consent prior to undergoing MRI when adhering to manufacturer conditions.

### Device interrogation and programming

Standard MRI safety protocols should always be followed for all patients, with the additional steps outlined below for scanning *MR Conditional* CIEDs (apart from implantable cardiac monitors). A full CIED interrogation should be performed to identify any device problems or higher-risk features that have not been identified previously or may lead to MR conditions not being met. This device interrogation can then be used for baseline comparison in the event of device parameter abnormalities being detected post-MRI.

The CIED should then be programmed, ideally just before the MRI scan. CIED interrogation and programming should be performed by an appropriately trained cardiac physiologist (or cardiologist), and all prescan parameters should be documented to enable comparison for changes postscan.

The programming changes for MRI should adhere to manufacturer instructions by using an MRI mode. Within MRI mode, there are generally programming choices to enable patients both with and without stable underlying rhythm to be scanned safely. Typically, pacing should be programmed off with OVO or ODO

modes if there is an adequate underlying rhythm; or asynchronous pacing with VOO or DOO modes if there is significant bradycardia or an unstable underlying rhythm.<sup>1 2 35</sup> Decision-making regarding the CIED mode to be programmed for MRI will be led by the cardiac physiologist or cardiologist. If there is doubt about the risk of competitive rhythm, it may be necessary to observe the rhythm for a few minutes after programming and discuss the programming options with a cardiologist. For all ICDs, antitachycardia therapies will be disabled in all MR modes. If CIED programming occurs in a different department to the MRI scan, the patient may require monitoring during transfer depending on the risk of the programmed settings.

For some older implantable cardiac monitors, data should be downloaded prior to the MR scan as this may be corrupted following exposure to the strong magnetic fields. Where feasible, a record of the download should be made in the patient's notes to alert radiography staff that this has been performed.

### Scan protocol

The relevant MR healthcare professional should check that an appropriate comprehensive protocol is provided prescan both to minimise the risk of requiring repeat scanning and ensure that MRI scanning conditions are met. When scheduling and coordinating scans, all appropriate disciplines should be available at the time of the scan.

### During scan

#### Monitoring during the scan

As with any patient, the risk of adverse events within the MRI scanner remains for the duration of the exam, and patients with CIEDs should be managed no differently. In the event of cardiac arrest or anaphylaxis due to contrast administration, the patient should be treated according to standard MRI department operating procedures. All staff must be familiar with evacuation protocols from an *MR Environment*.

Verbal communication with the patient is strongly recommended where possible. There is a low correlation of patient-reported symptoms and objective end points of lead parameter changes or MRI parameters associated with theoretically increased risk.<sup>35</sup> Nevertheless, because the nature or timing of an event is unpredictable, communication for new symptoms or unresponsiveness may provide an early sign of deterioration.

Pulse waveform monitoring (plethysmography) is the recommended minimum method of monitoring heart rhythm, and is generally resistant to artefact from electromagnetic interference (EMI) during the scan. This can either be obtained via pulse plethysmography sensors built into MR scanners or via dedicated *MR Conditional* monitors. ECG or vector cardiography (VCG) monitoring is an important additional method that is recommended where feasible, and may be required to meet the requirements of some *MR Conditional* CIEDs. While some manufacturer guidance for *MR Conditional* CIEDs state only ECG monitoring is required, it may be difficult to monitor acute rhythm changes if there is artefact on the ECG signal from the gradient fields. Blood pressure monitoring is not recommended as an alternative to ECG and pulse waveform monitoring. While this was performed in the MagnaSafe registry, cycling of a blood pressure cuff is likely to introduce delays to definitive assessment of an acute rhythm disturbance.

There should be a nominated physician responsible for the safety aspects available within the hospital at the time of scanning who is aware that a CIED MR is being performed. This may be a radiologist or cardiologist. There should also be both

someone able to reprogramme the CIED and personnel who are adult cardiac life support (ACLS) trained available in the hospital at the time of the scan. For inpatient National Health Service facilities, the resuscitation team can generally provide the ACLS trained staff, but for other models of care specific arrangements will need to be made.

One staff member must be present who is able to monitor and detect a change in the patient's rhythm from the available monitoring equipment. This could be a trained radiographer, trained nurse, cardiac physiologist or physician, provided they are capable of recognising significant changes in the heart rate/rhythm. The choice of staff member should be made locally, based on the knowledge and experience of the staff.

We recommend that staff with basic life support accreditation or above are present within the MRI unit for the duration of the scan. We recognise that other guidelines state a clinician with ACLS training is present with the patient from initial device programming to reprogramming after the MRI scan.<sup>1</sup> These recommendations go beyond manufacturers' guidelines, and some CIEDs permit auto-mode switching up to 48 hours prior to the MRI scan, and for up to 48 hours afterwards.<sup>36</sup> Resources for CIED implantation or revision do not need to be available on the same site as the MRI facilities when scanning *MR Conditional* CIEDs.

### Scan acquisition

At the time of writing, the most conservative conditions for *MR Conditional* CIEDs typically only require the MR Operator to ensure the MR scanner is in normal operating mode with regard to the specific absorption rate (SAR). There may be additional conditions including patient positioning, exclusion zones and field strength limits, although these can change as further manufacturer testing is performed. MR conditions often preclude the use of local transmit-receive coils directly over the CIED, such as a transmit-receive 31P-coil, but typically such coils are only used in research studies. The scan protocol should be decided prior to the patient arrival to include only the sequences required for a fully diagnostic examination. Metal artefact reduction strategies may be required to obtain diagnostic imaging, but are generally not required, especially for non-thoracic scans.<sup>37–39</sup> To avoid patient recall and associated logistical issues, each radiology department should have a mechanism in place to check that the images acquired are diagnostic and sufficient for reporting, prior to scan completion.

### Terminating a scan

In the event of a suspected arrhythmia, it is the responsibility of the attending MRI staff to evacuate the patient from the *MR Environment* as quickly as possible. Evacuation from the scanner is typically coordinated by the radiographer. One member of the team will administer basic life support while the cardiologist/cardiologist interrogates the CIED or applies an external magnet to the generator. In the case of a bradycardia, the CIED can be programmed to pace asynchronously. In ventricular arrhythmias, the ICD can be programmed to deliver appropriate therapy. If there are any delays in CIED interrogation, ACLS protocol should be followed with the attendance of the cardiac arrest team.

### Postscan

The CIED should be re-interrogated and programmed back to its original settings. If there is any significant change in a parameter, this should be communicated to the patient and

**Table 2** Summary of roles and responsibilities of each team member regarding MRI scanning of patients with *MR Conditional* CIEDs

Person	Key roles and responsibilities for MRI with <i>MR Conditional</i> CIEDs
Patient	<ul style="list-style-type: none"> <li>▶ Ensure that the referrer and MRI departments are aware of presence of CIED prior to attending the scan.</li> <li>▶ Facilitate provision of CIED information to radiology departments prior to scanning.</li> </ul>
Cardiologist/Cardiac physiologist	<ul style="list-style-type: none"> <li>▶ Ensure that <i>MR Conditional</i> CIED systems are implanted wherever possible (including consideration of pre-existing leads during generator exchange).</li> <li>▶ Include discussion regarding MR labelling in formal consent pre-implant (essential for <i>MR Unlabelled</i> CIED or <i>MR Conditional</i> CIED outside specified MR conditions).</li> <li>▶ Patient education postimplant regarding MRI and who to contact if difficulties accessing MRI.</li> <li>▶ Document <i>MR Conditional</i> labelling of the system and each component in medical records and on CIED identification card.</li> <li>▶ Communication with radiology department regarding CIED details, MR labelling and any potential exclusions (eg, abandoned leads), if requested during scan booking.</li> <li>▶ Device interrogation and reprogramming pre-MRI and post-MRI.</li> <li>▶ Rhythm monitoring during scan (dependent on local protocols).</li> </ul>
Referrer	<ul style="list-style-type: none"> <li>▶ Decision to perform and refer for MRI based on the same factors as for patients without CIEDs.</li> <li>▶ Identify the presence of CIED on MR request.</li> <li>▶ Provide clinical indication for the scan to enable appropriate protocol.</li> <li>▶ Liaise with cardiologist and patient to provide CIED details to MRI centre.</li> </ul>
<i>MR Clinicians</i> (radiologist/imaging cardiologist)	<ul style="list-style-type: none"> <li>▶ Establish process for accepting MRI referrals and identifying the appropriate radiologists who will check and report scans.</li> <li>▶ Review referral and prescribe scan protocol within MR conditions.</li> <li>▶ Patient safety during the MR examination.</li> <li>▶ May be required to check completeness of scan information before patient leaves the MR scanner.</li> </ul>
<i>MR Responsible Person</i>	<ul style="list-style-type: none"> <li>▶ Ensure appropriate local MR safety policies are in place.</li> <li>▶ Ensure MR staff have appropriate MR safety training.</li> <li>▶ Advise on the procurement and assessment of <i>MR Conditional</i> monitoring equipment.</li> </ul>
<i>MR Safety Expert</i>	<ul style="list-style-type: none"> <li>▶ Provision of MR safety advice.</li> <li>▶ Advise on MR sequence optimisation to meet scanning conditions and reduce artefact.</li> <li>▶ May assist in MR safety training.</li> <li>▶ Advise on the procurement and assessment of <i>MR Conditional</i> monitoring equipment.</li> </ul>
MR Radiographer	<ul style="list-style-type: none"> <li>▶ Check all local MR safety processes are followed for each patient.</li> <li>▶ Ensure MR conditions are adhered to.</li> <li>▶ Patient safety and communication before, during and after the MR examination.</li> <li>▶ Ensures that the scan is complete and will not require recall, or seeks advice where needed.</li> <li>▶ Rhythm monitoring during scan (dependent on local protocols).</li> </ul>
Most units will have several staff members trained to fulfil each of the roles, and there may be overlap in responsibilities. CIED, cardiac implantable electronic device.	

ongoing follow-up should be arranged. It is left to the discretion of the cardiology team what constitutes a change that is not due to physiological fluctuation and measurement imprecision, and the timing of follow-up.<sup>35 40</sup> All measurements should be documented.

It is also recognised that the beeper alarm function of CIEDs from several manufacturers may be permanently disabled after MRI (even in *MR Conditional* CIEDs), necessitating home monitoring or more frequent CIED follow-up particularly where devices are already under a manufacturer advisory. Cardiac physiologists should alert the patient and their usual device clinic if they feel that this loss of beeper function should lead to a change in routine follow-up protocols.

Table 2 summarises the key roles and responsibilities associated with individuals and departments regarding *MR Conditional* CIEDs. In some organisations, some of these responsibilities may be assigned to other individuals.

### Recommendations: infrastructure requirements

#### Personnel

Providing MRI for patients with CIEDs means that specialties need to work beyond traditional silos of practice. Staff involved in this service will include team members from the Departments of Radiology and Cardiology. Typically, this constitutes radiologists, MR radiographers, *MR Safety Experts*, cardiologists, cardiac physiologists. Other administrative and clerical

staff also play a vital role, and their additional time and need for training must also be recognised. It may be useful to accelerate expertise and communication using named individuals, at least when establishing services. Institutions should have clearly defined protocols in place to ensure that these general recommendations account for the local environment and variations in service requirement. Local protocols must provide details regarding escalation procedures in the event of complications, with procedures for instituting basic and advanced life support where needed.

#### MR scanner

Up-to-date information regarding the MR conditions for scanning the CIED should be followed. MR conditions for *MR Conditional* CIEDs are relatively simple to meet on all current clinical MR systems. All devices allow for scanning within cylindrical bore 1.5 T systems and often 3 T scanning is accepted within the CIED manufacturer conditions. Other variations in MR scanner hardware/software (eg, maximum gradient amplitude) typically do not present a limitation for scanning patients with *MR Conditional* CIEDs. Advice from the *MR Safety Experts* can be sought if required.



### Monitoring equipment

Patients must be monitored throughout the scans with a minimum of pulse oximetry waveform monitoring and/or ECG. The provision of *MR Conditional* monitoring equipment may require investment for units starting scanning patients with CIED, and sites procuring new *MR Conditional* monitoring equipment are encouraged to seek assurances that the equipment will provide robust data during MRI. Notably, most ECG monitoring systems—even many that are *MR Conditional*, are susceptible to significant signal artefact during some MR sequences. We therefore strongly recommend that continuous pulse oximetry waveform monitoring (which is generally resistant to artefact from gradient fields during image acquisition) is performed in all patients.

Although not approved for diagnostic purposes, MR scanners have their own monitoring systems with ECG and pulse oximetry waveform assessment for gating during image acquisition, and these may be sufficient for detecting changes in rhythm (rather than QRS/ST segment changes) needed for the purposes of monitoring patients with CIEDs. This approach has precedent given that it is currently routinely used during stress perfusion cardiac MRI, although this remains an off-label use of the scanner.

### CIED programming units

All CIEDs except implantable cardiac monitors require programming before and after MRI. This requires the availability of a pacing system analyser (PSA), which is a portable unit that is able to interrogate and programme CIEDs. These are specific to the manufacturer and are available in all cardiac physiology departments or via manufacturers, but all are *MR Unsafe*. For patients with *MR Conditional* CIEDs, there should be a PSA available within the hospital at the time of scanning, however this does not need to be physically in the MRI department. Patients can be programmed within the cardiology or pacing department prior to arriving at the MRI unit.

### Resuscitation equipment

A resuscitation trolley should be available within the MRI department. This should contain an external defibrillator. A manual defibrillator with the ability to transcutaneous pace may not be available in all MRI units, although is recommended when scanning patients with ICDs or those who are pacing dependent. For units scanning patients with *MR Unlabelled* devices, a manual external defibrillator with transcutaneous pacing capability must be available. Requirements should be discussed locally with the cardiology team. An external CIED magnet (available from CIED clinics) that can be applied to the CIED should also be available to enable reprogramming of the CIED to a default setting in an emergency situation.

## CONSIDERATIONS AND RISKS OF DEVICE REPROGRAMMING IN PATIENTS WITH *MR CONDITIONAL* CIEDS

For all patients with cardiac pacemakers, defibrillators or cardiac resynchronisation therapy devices, the CIED must be interrogated and reprogrammed into an 'MRI mode' prior to undergoing the scan in order to minimise the risk of complications from inappropriate pacemaker/defibrillator activation or inhibition of pacing. There are some additional considerations related to choice of CIED programming modes that should be adhered to in specific devices/patients.

### Potential risks associated with changing CIED programming mode

Prior to scanning, it is important to ascertain the underlying heart rhythm and pacing requirement for patients with CIEDs in situ, even for patients with *MR Conditional* devices. Pacemaker-dependent patients have no underlying heart rhythm, or their intrinsic heart rate is sufficiently slow to cause symptoms and make the patient haemodynamically unstable, should pacing not be delivered.

Without appropriate device reprogramming, in pacemaker-dependent patients there is a risk that the device interprets EMI from the MRI scan as spontaneous myocardial activity (oversensing true cardiac electrical activity) and in response inhibits pacemaker function. To minimise this risk of oversensing, devices should therefore be programmed either to pace continuously (asynchronously) for patients with high pacing demand or those who are pacemaker dependent with no underlying rhythm, or alternatively with pacing programmed off for the duration of the scan where the patients' underlying heart rhythm is stable.

For patients with an acceptable stable intrinsic heart rate and rhythm, there is a risk of arrhythmia either through unintentional pacemaker activation if the device is not reprogrammed into MRI mode, or with intentional programming to an asynchronous pacing mode that competes with the patient's intrinsic rhythm. Pacing at the same time as intrinsic cardiac repolarisation risks ventricular arrhythmia (termed 'R on T'), although asynchronous pacing is routinely performed during pacemaker lead threshold checks, with only extremely rare published cases of ventricular fibrillation precipitated by this (quoted risk <0.001%).<sup>41</sup>

### Potential risks to patients with defibrillators

Defibrillator antitachycardia and shock therapies need to be programmed off for the duration of the MRI scan for all patients with ICDs and cardiac resynchronisation therapy-defibrillators (CRTDs). If defibrillator antitachycardia therapies are left activated for the MRI scan, this may result in inappropriate therapies or device malfunction. EMI is likely to be interpreted by the CIED as ventricular tachycardia leading to attempted delivery of therapy in the form of antitachycardia pacing or shocks. For this reason, all ICD therapies are automatically programmed off in all MR modes in all *MR Conditional* devices.

This does however mean that if a patient were to develop a ventricular arrhythmia while in the scanner, they would need to be evacuated and treated outside the *MR Environment* with external defibrillation or by reactivating the defibrillation therapies via reprogramming. ICD implantation is performed in patients at increased risk of ventricular arrhythmias—either for secondary prevention (in survivors of cardiac arrest or frequent ventricular arrhythmias) or primary prevention for patients with underlying cardiac conditions predisposing them to high risk of arrhythmias. For a 1 hour MRI scan in a patient with a standard primary prevention indication for ICD implantation, the risk of ventricular arrhythmia during an MRI scan is approximately ~0.0005%.<sup>42 43</sup> Given this extremely low risk, patients with *MR Conditional* ICDs are treated as having a similar risk profile to patients with *MR Conditional* permanent pacemakers. If a recent clinically significant ventricular arrhythmia is detected prior to imaging, the cardiologist should be consulted to provide an opinion as to the risk of proceeding with the scan and programming the tachycardia therapies off.



### Potential risks to patients with cardiac resynchronisation therapy pacemakers

Cardiac resynchronisation therapy-pacemakers provide biventricular pacing for patients with heart failure to maintain synchronous myocardial contraction, improve cardiac output and alleviate symptoms. Currently, 'MRI-mode' for most cardiac resynchronisation therapy-pacemakers permits only right ventricular (rather than biventricular) pacing, resulting in temporary loss of biventricular pacing. Although cardiac output will fall to a small extent while the device is programmed without biventricular pacing, the short duration of the scan means that the clinical risk of cardiac decompensation is negligible in haemodynamically stable patients undergoing MRI.

### Potential risks to patients with implantable cardiac monitors

Implantable cardiac monitors include implantable loop recorders and pulmonary artery pressure monitors. There has been no reported harm to patients with implantable cardiac monitors undergoing MRI according to manufacturer stated conditions, and all modern devices are *MR Conditional*. For some implantable cardiac monitors, it is recommended that data are downloaded from the devices prior to undergoing MRI, but this is not needed for newer devices. Importantly, interrogation and reprogramming pre-MRI and monitoring during scans are not needed for patients with implantable cardiac monitors.

### MR UNLABELLED CIED SYSTEMS AND MR CONDITIONAL CIEDS NOT FULFILLING SPECIFIED CONDITIONS

Alongside CIED systems without MR safety labelling, there are many different scenarios where one or more of the specified conditions for an *MR Conditional* CIED are not met, each of which are associated with different levels of MR safety risk. Importantly, although the risks may be greater than when scanning patients with *MR Conditional* CIEDs where all the conditions are met, for many scenarios the lack of MR safety labelling reflect the limitations of testing, namely a lack of sufficient evidence to demonstrate that the device is safe under certain conditions, rather than any evidence demonstrating that it is unsafe. If the benefit to the patient outweighs these risks and there are no alternative modalities able to answer the clinical question, then the scan should be carried out, provided appropriate steps (as outlined below) are performed and documented. Risks should be mitigated where possible and the clinical decision should take into account the risks regarding onward clinical management if a decision is made not to perform the MRI

(eg, both the risks of invasive biopsy and that of diagnostic uncertainty from incorrect or incomplete diagnosis).

For the many patients who find themselves in these scenarios, access to MRI currently can be particularly challenging.<sup>15</sup> Given the availability of increasing safety data and the high clinical need for individual patients, some centres with good collaborative working between cardiology and radiology departments may wish to provide MRI scans for these patients. We aim to encourage a network of centres to cater for regional demand across the UK. This will permit centres to gain from efficiencies of scale, and centralisation of expertise.

The recommendations below are consistent with MHRA guidelines for scanning patients with implants where MRI may be contraindicated.<sup>23</sup> figure 3 describes the suggested workflow when considering an MR request for a patient with an *MR Unlabelled* device or an *MR Conditional* CIED where one or more of the conditions cannot be met.

### Risk assessment and risk-benefit analysis


A risk assessment and risk-benefit analysis should be undertaken with involvement from a combination of the radiologist, cardiologist, referrer, *MR Operator*, *MR Responsible Person* and the *MR Safety Expert*.

The following points should be confirmed and documented prior to a decision to scan:

1. The MRI scan is likely to change patient management.
2. There is no appropriate alternative modality to answer the clinical question.
3. The potential benefit outweighs the risk of the MRI scan.

Risk-benefit analyses should incorporate individual patient, device and scan-related factors with specific scenarios as outlined below. The risks associated with a particular scenario can vary significantly depending on the individual circumstances, and consensus opinion is provided where clinical evidence is limited.

Metallic artefact from the CIED generator can provide additional challenges for anatomical regions that lie close to the device. This is a particular problem for cardiac MRI studies, however diagnostic imaging is generally feasible with published strategies.<sup>39 44</sup> Artefact is generally more frequently encountered with ICD and CRTD devices, and *MR Conditional* labelling does not attenuate this problem or guarantee diagnostic image quality. Similarly artefact tends to increase at higher field strengths, hence imaging at 1.5 T rather than 3 T is recommended where

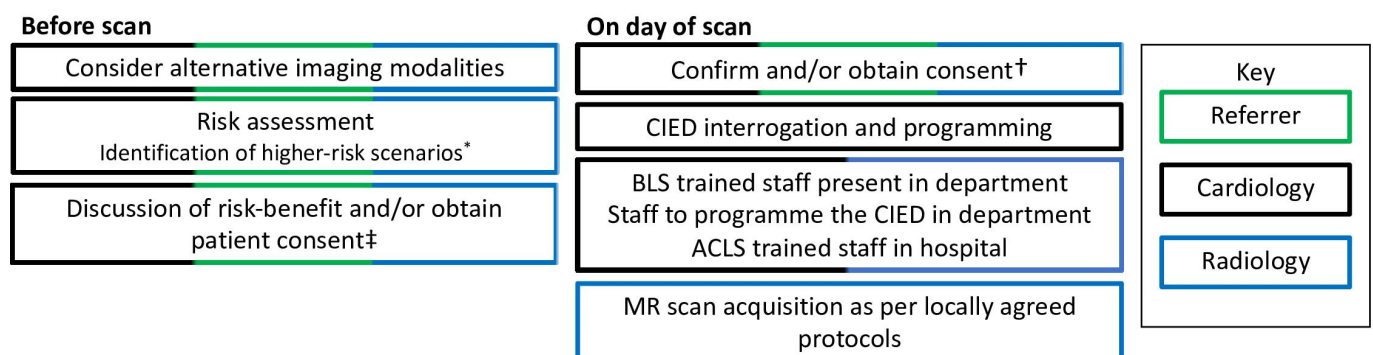

**Figure 3** Suggested additional steps required prior to performing MRI for patients with CIEDs that are MR Unlabelled or do not fulfil MR conditions.

\*Higher-risk scenarios include the presence of fractured, epicardial, abandoned leads; recent implantation; battery at ERI; deactivated systems; lead parameters outside manufacturer recommendations, other implants present. †Appropriate person confirming consent decided as per local protocol.

‡Consent can take place on the day of the scan, according to local protocols. For the purposes of this document, Present=present at scanner side.

ACLS, adult cardiac life support; BLS, basic life support; CIED, cardiac implantable electronic device; ERI, elective replacement indicator.

## Guideline or consensus statement

**Table 3** Risk stratification of performing MRI in patients with MR Unlabelled CIEDs or MR Conditional CIEDs outside specified conditions

Risk category	Scenario	Clinical risk	Level of evidence
Lowest	MR Conditional CIEDs* (meeting all conditions)	MR safety tested by device manufacturer. <sup>25–27</sup>	A
	MR Conditional CIEDs with additional MR Conditional device implanted (eg, coronary stent)	No clinical evidence of increased risk.	A
Lower	Recent implants (<6 weeks)	No clinical evidence of increased risk. <sup>49</sup>	C
	Unmet condition due to presence of additional implanted device	No clinical evidence of increased risk.	C
	Temporary surgical epicardial pacing wires (with no external component)	No clinical evidence of increased risk. <sup>66–69</sup>	C
	Unmet patient position exclusion zone, MR Conditional CIEDs	No clinical evidence of increased risk. <sup>78 79</sup>	B
	Scanning beyond SAR restrictions with MR Conditional CIEDs when required for diagnostic imaging	No clinical evidence signal of increased risk. <sup>75 76</sup>	C
	3 T MRI field strength, MR Conditional CIEDs labelled as MR Conditional at 1.5 T only	No clinical evidence of increased risk. <sup>80 81</sup>	C
	MR Unlabelled pin plug with MR Conditional generator and leads	No clinical evidence of increased risk.	C
	'Mismatched' CIEDs with MR Conditional generators, MR Unlabelled leads	No clinical evidence of increased risk. <sup>34</sup>	C
Intermediate	Inactive, battery-depleted CIEDs	Potential signal of increased risk from one case. <sup>34</sup>	C
	Generators implanted outside the pectoral region	Higher risk.	C
	MR Unlabelled generators, any lead MR labelling (non-pacing dependent)	Risk of device failure.	B
	MR Unlabelled generators, any lead MR labelling (pacing dependent)	Higher risk of asystole in the event of device failure.	B
	Abandoned leads (capped or not)	Lead/Tissue heating in experimental studies, no reported clinical complications. <sup>49 57 59 60 63–65</sup>	C
Higher	Stable abnormal lead parameter	Higher risk, mitigated by investigating cause and appropriate device programming.	C
	CIED component advisory warning	Higher risk depending on cause.	C
	Permanent epicardial leads	No clinical evidence of increased risk. <sup>49 67</sup>	C
	Fractured leads	Higher risk.	C
	Battery at elective replacement interval	Higher risk.	C
	Temporary systems with externalised generator	No clinical evidence of increased risk. <sup>49 71</sup>	C
	Pre-2000 market release generators	Increased risk of electrical reset. <sup>45</sup>	B
Avoid	Pacing-dependent patients with ICDs devices where not possible to maintain asynchronous pacing (VOO/DOO)	Unavoidable extremely high risk—do not scan.	C
	Scanning any active CIED without reprogramming prior to scan to minimise risk	Risk highly possible—do not scan.	C

Risk is represented as categories, but within each category the scenarios should be considered as a scale from low to progressively higher risk (top to bottom of the table - shaded from green to red respectively, with black shaded categories representing scenarios where MRI should be avoided). Overall risk is a combined consideration of the likelihood of an event and the clinical outcome if that event occurred. Scenarios considered higher risk either have a paucity of trial data describing clinical safety end points, or are based on expert consensus despite small series reporting safe MRI scanning. Level of evidence A=data derived from multiple randomised controlled trials, or meta-analyses; level of evidence B=data derived from a single randomised clinical trial or large non-randomised study; level of evidence C=consensus of expert opinion based on clinical experience or case series.

\*Includes both permanent pacemakers. ICDs and implantable cardiac monitors.

CIED, cardiac implantable electronic device; ICD, implantable cardioverter-defibrillator; NA, not applicable; SAR, specific absorption rate.

available. For most MRI scans, artefact should not be a major factor in decision-making.

### Risks associated with MRI in specific scenarios

A summary of the risks associated with the scenarios discussed below is provided in [table 3](#).

#### MR Unlabelled CIEDs

The risks associated with any MR Unlabelled active implantable medical devices include mechanical forces (attraction, torque, vibration), heating, unintended stimulation and device malfunction. CIED generators also may contain components sensitive to the magnetic field including reed switches, which may change position in an MR field leading to alterations in programming mode. Software corruption can occur due to electromechanical interference (EMI) known as 'power-on reset', causing the CIED to revert to a back-up mode of programming. Before the introduction of MR Conditional CIEDs, developments associated with CIED design had already gone some way to reducing these

risks. Even for pacemakers with market release before 2002, the attractive force and torque experienced due to the magnetic fields associated with a 1.5 T MRI scanner were shown to be lower than the confining forces of surrounding tissue and hence low enough to present no safety risk.<sup>45</sup> Technical developments associated with CIED generators generally (particularly ICDs) mean those implanted after approximately the year 2000 have reduced risk of heating, malfunction and electrical reset during MRI both in testing and from clinical data.<sup>45</sup> Implantation dates before 2000 and 2002 were used as exclusion criteria for the two largest studies of MR scanning of patients with MR Unlabelled CIEDs.<sup>35 46</sup>

There is a significant and growing body of clinical evidence to support MR scanning of patients with MR Unlabelled CIEDs under strictly controlled conditions, similar to MR Conditional CIEDs. Three recent registries totalling 2859 patients with MR Unlabelled CIEDs undergoing MRI have reported no deaths or life-threatening arrhythmia.<sup>35 47–49</sup> Patients were however excluded if they were pacing-dependent and had an ICD without

asynchronous pacing capability. Fifteen of 2859 patients with CIED underwent partial or complete software resets ('power-on resets'), many of which were transient or could be programmed around. One patient required a generator change as the patient did not undergo CIED reprogramming prior to an MRI scan.<sup>50</sup> A recent meta-analysis of 5625 patients undergoing 7196 MRI scans similarly reported no deaths and no lead complications. There were 1.4% of cases with power-on resets, but none in devices that were released to market after 2005.<sup>51</sup> Another meta-analysis of 5099 patients (overlapping with these studies) included one report of inappropriate ICD shock delivery (n=1).<sup>49</sup>

Although the number of patients included in published data is relatively large, it should be noted that the number of cases with a particular combination of generator and leads may be very small. To satisfy regulatory bodies that a particular combination can be labelled as *MR Conditional* typically involves modelling millions of such potential exposure conditions.<sup>52-53</sup>

These data highlight the low risk of MRI scanning, provided that strict protocols are followed, although quantifying personalised absolute risk is complex because of individual differences in pacing and device component factors. In summary, there is a significant body of clinical evidence to support MR scanning of patients with *MR Unlabelled* CIEDs implanted after 2002 under strictly controlled protocols similar to *MR Conditional* CIEDs, when clinically indicated. Patients with *MR Unlabelled* CIED generators with market release prior to this time appear to have a slightly greater risk from MRI including around 1.5% risk of electrical reset alongside other unpredictable complications,<sup>54</sup> and therefore this should be considered in the risk-benefit assessment prior to scanning—particularly for patients with high pacing requirements. MRI is not recommended for patients who are pacing dependent and have an ICD without asynchronous pacing capability.

#### *MR Conditional* CIEDs with 'mismatched' CIED components

CIED generators and leads are manufactured and sold as separate components, meaning that implantation of a fully *MR Conditional* system requires operator selection of appropriate individual components. 'Mismatched' CIED systems are those with either only some components which are *MR Conditional*, or with fully *MR Conditional* components but produced by different manufacturers, and so will not have been formally tested in combination.

Manufacturers have performed formal MR safety testing of older leads that were previously *MR Unlabelled*, and have frequently shown them to be sufficiently safe when combined with *MR Conditional* generators from the same manufacturer to satisfy retrospectively relabelling the leads as *MR Conditional*. Comprehensive formal testing of every possible generator-lead combination (especially between different manufacturer components) is neither feasible nor appropriate meaning that this issue is unlikely to be eliminated.<sup>55</sup>

Patients with *MR Unlabelled* generators and leads are likely to require generator exchange before the battery reaches end of life (EOL) (5–10 years for the majority of CIEDs), and currently the majority will have a new *MR Conditional* generator implanted and connected to the original *MR Unlabelled* leads. Recent multi-centre data have found no increased incidence of adverse effects of MRI with *MR Unlabelled* leads as compared with *MR Conditional* leads—both in terms of safety events and changes to lead parameters.<sup>34</sup> These data (in combination with that from other studies with *MR Unlabelled* devices) suggest that the clinical

risk of MRI in patients with *MR Conditional* generators is not increased by having *MR Unlabelled* leads connected. To facilitate equitable access to MRI, many centres currently scanning only *MR Conditional* CIEDs may choose to regard MR scanning for 'mismatched' devices with *MR Conditional* generators as a lower risk scenario so standard protocols for *MR Conditional* CIED systems are used, provided that there are no other high-risk features.

#### *MR Unlabelled* pin plugs

There are clinical scenarios when a patient does not require implantation of all of the leads that can be accommodated by a particular generator. In these situations, a pin plug is inserted into the generator to fill the port at the time of implantation. This serves the purpose of blanking off the port to ensure biological tissue does not enter and should not affect the electrical configuration of the device. Plug attachment types follow an industry standard and hence this compatibility between manufacturers can lead to implanted systems with a mismatched port plug. To form a complete *MR Conditional* pacing system, specific pin plug models have been tested by manufacturers in combination with other components. Consequently, use of pin plugs from a different manufacturer may invalidate an *MR Conditional* of the *MR Conditional* generator. However, these are considered a very low-risk scenario with no reported adverse effects, provided there are no other high-risk features.

#### Non-standard lead implants and additional leads

CIED leads are generally implanted permanently via venous access to the heart. However, leads can also be implanted with the intention of remaining temporarily, or may be implanted by differing access including surgically implanted leads that are attached to the epicardium of the heart. CIED leads can also malfunction, fracture or dislodge over the many years that they are implanted. To ensure ongoing device function and given the risks of lead extraction, additional leads may be implanted and the non-functioning lead removed from the generator but left in situ. A functional or non-functional lead that is left in place and is not connected to a CIED is termed an abandoned lead.<sup>29</sup> If there is any doubt regarding the implanted hardware and associated risks, patients should undergo chest X-ray and there should be further discussion with a cardiologist with appropriate experience of CIEDs and MRI. The risks associated with MR scanning of patients with abandoned or fractured leads include the potential for induced voltages in the leads from the RF field (strongest in amplitude within the transmit RF coil) or the imaging gradients (strongest in amplitude around 30 cm away from magnet isocentre), causing lead heating (most likely at the lead tip) and/or direct stimulation of the cardiac muscle. One method to mitigate the risk from the RF field is through the use of transmit receive (T/R) coils. These ensure the majority of the RF energy is only imparted to a specific anatomical region within or covered by the coil and therefore, if positioned appropriately, can minimise the RF exposure to the CIED.

#### *Abandoned or fractured transvenous permanent pacemaker or defibrillator leads*

Experimental evidence to help quantify the risks of abandoned or fractured transvenous permanent pacemaker and defibrillator leads is mixed with some studies suggesting large temperature increases both in vitro<sup>56-57</sup> and in vivo animal studies.<sup>58</sup> More recently, in vitro testing at 1.5 T has demonstrated greater MRI-induced heating in abandoned leads (up to 29.9°C) compared



with pacemaker-connected leads (up to 11.6°C),<sup>59 60</sup> although temperature rises were strongly dependent on lead length and were generally higher when the abandoned leads were capped. A recent study performing electromagnetic simulations of numerical models with fractions of retained endovascular leads positioned at different imaging landmarks found temperature rises during a 10 min scan at all imaging landmarks remained <3°C at 1.5 T and <6°C at 3 T.<sup>61</sup> Another numerical study at 1.5 T found the deposited power at the lead tip for fractured leads (increased up to a 16 times compared with non-fractured leads) was dependent on type of conductor break, and the design and location of the lead.<sup>62</sup> Importantly, in vitro measurements demonstrating temperature rises up to 7°C around the lead tip were found to be insignificant when repeated in vivo, where there is a significant cooling effect from adjacent blood flow.<sup>45</sup> Additionally, any localised RF heating is expected to be significantly less when the device/leads are located outside of the transmit RF coil.

Despite these experimental data, clinical reports have not described clinical or electrical evidence of CIED dysfunction, arrhythmia or pain as a result of MRI in patients with abandoned leads. Clinical studies have reported outcomes on patients with either abandoned leads alone, or in the presence of an additional functioning CIED, and so both scenarios are included in this recommendation. Currently, the largest clinical dataset available included 200 scans in 139 patients with CIEDs and abandoned leads, with no clinical complications, similar to other case series.<sup>63–65</sup> Given the mixed evidence from simulations and preclinical work together with the relatively small clinical evidence of scanning patients with abandoned leads, this scenario is considered intermediate risk.

#### *Permanent epicardial pacing leads*

To our knowledge, currently no permanent (surgically implanted) epicardial pacing leads are labelled *MR Conditional*, due to the risk of RF heating at the lead tips, which has been demonstrated ex vivo.<sup>49 59 60 63 64</sup> For these reasons, patients with surgically implanted permanent epicardial leads were excluded from large registries performing MRI scanning of patients with *MR Unlabelled* CIEDs. However, to our knowledge no clinical adverse events have been reported from scanning CIEDs with the presence of epicardial leads, and many units currently scan such devices using the additional precautions recommended for patients with *MR Unlabelled* devices.<sup>59 64</sup> All patients with permanent epicardial pacing leads should be considered as a high-risk MRI scenario even if other components are *MR Conditional*.

#### *Postoperative epicardial pacing wires*

Temporary epicardial wires placed at the time of cardiac surgery are different to surgically implanted epicardial leads as part of a permanent CIED system. Where possible, temporary epicardial wires will be removed postoperatively. Alternatively, these may be cut at the skin after surgery leaving a short length of wire implanted chronically. A study including 51 patients who underwent MRI at 1.0 or 1.5 T with temporary epicardial pacing wires cut short at the skin found no reports of clinical events or symptoms suggesting arrhythmia or other cardiac dysfunction.<sup>66–69</sup> Patients with postoperative epicardial pacing wires are fairly common, and given the low risk of complications, this should not be considered a contraindication to MRI scanning.

#### *Temporary CIED system with externalised generator ('temporary-permanent CIED')*

These are systems that are implanted temporarily (generally while patients are treated for systemic infection or where recovery of intrinsic electrical conduction is expected) and consist of a generator fixed to the external chest wall with a transvenous (usually active fixation) lead attached and implanted internally. There are few reports of such devices undergoing MRI,<sup>70 71</sup> and although clinical complications have not been reported, these should be considered for a high risk MRI, and this should only be performed when MRI is considered essential to the patient's clinical pathway.

#### *Recent device implantation (within minimum duration specified by MR conditions)*

Most *MR Conditional* CIEDs include the condition that scanning be performed a minimum of 6 weeks following implantation. This is often described as a period to allow fibrosis at the lead-myocardial interface. However, the theoretical risk of lead displacement is minimal given the lack of ferromagnetic components within *MR Conditional* leads, and hence negligible force on the lead tip.<sup>1</sup> Clinical scans performed within 6 weeks have not been associated with complications,<sup>72</sup> and similarly, no correlation between changes in lead performance (sensing, pacing threshold or impedance) and time from implantation was observed in cases from the MagnaSafe registry that included 17 cases in which MRI was performed within 30 days of implant, and 5 cases in which MRI was performed within 7 days of implantation.<sup>35</sup> Some manufacturers now provide flexibility in their MR conditions for scanning during this period if clinically necessary. It is therefore appropriate to perform an MR scan earlier than recommended if the scan indication is required, but patient positioning should avoid arm elevation above shoulder level for the first week postimplant as per standard postimplant care to reduce the risk of traction on the lead resulting in displacement.

#### *MR Conditional lead parameters that do not meet specified conditions*

Manufacturers of *MR Conditional* pacemakers and ICDs generally stipulate device parameter conditions that must be fulfilled prior to patients undergoing MRI, such as measured lead threshold and impedance values. For CIEDs with parameters outside of these values, the cause of any abnormal lead parameter should be investigated before MRI scanning by the cardiac physiologist and cardiologist. If there is evidence of lead fracture, this should be considered a high-risk MR scan. For patients with high lead thresholds or low sensing, if these abnormalities are stable and a sufficiently increased safety amplitude window can be programmed (at least twice threshold), systems should be treated as *MR Unlabelled* CIEDs and patients consented for relatively higher risk of complications.

#### *MR Conditional CIED generators with batteries close to depletion (elective replacement indicator or end of life)*

Lead outputs are significantly increased as default with MRI modes (generally to 5.0 V at 1.0 ms), and a drop in battery voltage has been observed for some *MR Unlabelled* CIEDs undergoing MRI. Together this means that manufacturers of *MR Conditional* devices stipulate a minimum battery voltage that generators should have prior to scanning in order to

fulfil conditions.<sup>73</sup> Excessive battery depletion during MRI of devices with prescan low voltages (elective replacement indicator (ERI) or EOL) may risk CIED malfunction or automatic pacing mode switch if the device reaches EOL. More recently, a review of 9 non-pacing-dependent patients with *MR Unlabelled* CIEDs who underwent 13 MRI within 3 months of the ERI reported electrical reset in 2 patients, although both occurred in pacemakers implanted pre-2005.<sup>74</sup> Generator change may be performed prior to MRI, where there is concern about the risks of mode switch or further battery depletion are high. Patients with batteries close to depletion should be considered to be at a higher risk of complications from MRI.

#### *Inactive, battery-depleted generators*

CIEDs may remain implanted with a depleted battery in patients without pacing requirements, where the clinical need for generator exchange is low. Although non-functioning, these devices are generally older and may therefore potentially still be at risk from MRI due to mechanical forces, heating, unintended stimulation. One case has been reported of tachycardia and chest pain on scan initiation within one large retrospective cohort (n=1148 MRI examinations).<sup>34</sup> Patients with inactive, depleted generators should therefore be considered as an intermediate risk for MRI.

*MR Conditional generators implanted outside of the pectoral region*  
CIEDs are generally implanted in the pectoral region but may (rarely) be implanted in other locations, including the abdomen. MR conditions generally stipulate that the generator should be implanted in the pectoral region, and the risks of MRI in patients with CIEDs in other anatomical locations are unlikely to have been tested by manufacturers. This means that such CIEDs, even if *MR Conditional*, should be considered as equivalent to *MR Unlabelled* devices, although the absolute clinical risk is not known. Given that patients may have abandoned hardware from previous implants in other locations, there may be other factors that increase the risk of MRI (eg, epicardial leads with an abdominal implant). Additionally, it is important to note any MRI exclusion zones will have been defined based on the specified device location.

#### *Unmet condition for additional implanted devices*

It is not possible for the manufacturer of an *MR Conditional* CIED to assess the potential interactions with all additional implanted devices. However, it is important to recognise that the risk of potential interactions between devices drops off significantly with separation distance, for example, ISO/TS 10974 only suggests an assessment for potential proximity enhancement from coupling between multiple electrodes is required when the separation distance is <2 cm. Although some *MR Conditional* CIEDs exclude the presence of any additional devices, others do not which suggests the risks associated with this scenario, particularly when the additional devices are well separated from the CIED, are low.

#### Scanner-related scenarios

##### *Unmet condition for MRI field strength or MRI scanner type*

The vast majority of clinical MRI scanners are closed bore cylindrical MRI systems operating at 1.5 T or 3 T. Although some *MR Conditional* CIEDs are labelled as *MR Conditional* at 1.5 T only, many now permit scanning at both 1.5 T and 3 T, although to our knowledge all are specified only for closed bore cylindrical

MRI systems. There is little evidence available for MRI scanning of patients with CIEDs at other field strengths or on other MRI scanner types, for example, open MRI systems. Importantly, it is unlikely that the different RF frequencies associated with other field strengths will have been assessed by the CIED manufacturer and therefore this scenario presents an unknown risk in terms of RF heating, device malfunction and unintended cardiac stimulation. For sites unable to meet the conditions of an *MR Conditional* CIED for MRI field strength and MRI scanner geometry, onward referral to a centre that can meet this condition is likely to be the most practical solution, although centres scanning *MR Unlabelled* CIEDs may choose to scan locally with the additional steps recommended for scanning patients with *MR Unlabelled* devices.

##### *Unmet condition for specific absorption rate (SAR)*

Many *MR Conditional* CIEDs require the MRI scanner to be restricted to the normal operating mode for SAR (whole body SAR limited to 2 W/kg). In general, this is achievable for many clinical MR sequences without significantly impacting image quality, but there may be occasions where there is a clinical need to operate in the first level controlled operating mode for SAR (whole body SAR limited to 4 W/kg). There is growing evidence to support safe MRI scanning of patients with CIEDs at these SAR levels. A study of 1464 patients with non-*MR Conditional* CIEDs who underwent 2028 MRI examinations without SAR restrictions found no evidence of an association between RF energy deposition, dB/dt, or scan duration and changes in device parameters.<sup>75 76</sup> Consequently, for scenarios where there is a clinical need to operate with SAR levels above the normal mode, the incremental risk appears to be relatively low.

##### *Unmet patient positioning exclusion zone condition or with thoracic isocentre*

Some *MR Conditional* CIEDs include an exclusion zone, to avoid positioning the device such that it is exposed to the highest levels of RF during the MRI scan. This may reflect conservative conditions that were incorporated into some clinical trials supporting regulatory approval of *MR Conditional* CIEDs, rather than confirmed evidence of risk. Indeed, many *MR Conditional* CIEDs no longer include this condition, suggesting the associated risks are low.

##### *Unmet patient decubitus condition*

Some *MR Conditional* CIEDs provide limited conditions for the patient decubitus during the MRI scan, for example, supine or prone. This is likely to reflect the limitations of MR safety testing and the simulations performed. To our knowledge, there have been no increased risks identified with different patient positioning and given several manufacturers of *MR Conditional* CIEDs do not state any such *MR Conditional*, the expected risk is low.

#### **Recommendations: additional protocol for *MR Unlabelled* CIEDs or *MR Conditional* CIEDs outside the specified conditions**

##### **Prescan preparation**

##### **Referral**

Patient and device factors that identify patients to be at higher risk for undergoing MR scanning should be identified before booking a scan (table 3). In addition to the particular risks posed in the *MR Environment* by CIEDs, it is important to note that CIEDs should be viewed as only part of the whole. Other MR

## Guideline or consensus statement

**Table 4** Suggested statements to use when describing risk during consent for patients with different non-MR Conditional CIEDs

Intervention	Recommended risk statement to discuss with the patient. <i>The MRI procedure, benefits and alternatives should also be discussed with the patient with the opportunity for them to have additional queries addressed by an appropriate clinician.</i>
<b>Intermediate-risk and higher-risk scenarios (formal written consent required)</b>	
<b>MRI with MR Unlabelled CIED (without additional higher-risk scenarios)</b>	You have been referred for a MRI scan. Your pacemaker/defibrillator has not been formally approved to undergo MRI scanning by the manufacturer. The decision to perform the MRI scan has been made after discussing the possible benefits, risks and alternatives with your referring doctor. Serious complications related to MRI occur in <1 in every 2000 patients (about 0.05%) with these devices overall. These include, but are not limited to: ► damage to the cardiac device; ► abnormal heart rhythms; ► excessive tissue heating. Emergency or urgent replacement of the cardiac device may be needed and will be performed if required.
<b>Additional intermediate-risk and higher-risk scenarios (formal written consent required)</b>	
<b>MRI with MR Unlabelled CIED generators implanted prior to 2005</b>	(in addition to above) Due to the age of your device, the risk may be slightly higher—with approximately a 2% risk of (generally temporary) programme changes to 'factory settings'.
<b>MRI with MR Unlabelled CIEDs implanted prior to 2000</b>	(in addition to above) There is less evidence for scanning patients with old devices that were implanted before the year 2000. We also know that the older technology used in these devices mean that they are more sensitive to MRI and therefore the risk is likely to be higher.
<b>MRI with abandoned lead(s)</b>	(in addition to above) Having a pacemaker or defibrillator lead which is not attached to a generator may result in heating at the end of the lead in your heart, which could theoretically cause tissue damage. To date, there have been no reported problems in patients being scanned with these leads, although the number of these patients is relatively small. We would ask that you inform staff immediately if you feel any discomfort.
<b>Lower-risk scenarios (formal written consent is not required, and the risks can be discussed verbally)</b>	
<b>MRI following recent CIED implantation (typically &lt;6 weeks post implant)</b>	Your cardiac device manufacturer recommends that you wait for a period of time after implantation (commonly 6 weeks) before having an MRI scan. There have been no problems reported in patients having scans earlier than this, however formal testing has not been performed to guarantee that this is safe.
<b>MRI in patients with 'mismatched' CIEDs with MR Conditional generators</b>	The generator in your device has been formally tested and approved for MRI, however the leads have not. Studies have found no increased risk of MRI with devices like yours when compared with MRI in patients with device systems that are fully approved. There is however a potential risk of undergoing an MRI, but this will be small and considerably lower than <1 in every 2000 patients, which is the overall risk of MRI in patients with devices not approved for MRI.
These statements should be used in addition to discussing the MRI procedure, potential benefits and alternatives. This list is intended for common scenarios, and not as an exhaustive list. 'Mismatched' CIEDs have MR Conditional generators and MR Unlabelled leads; or MR Conditional components from different manufacturers. *For lower risk scenarios, see <a href="#">table 3</a> . CIED, cardiac implantable electronic device.	

contraindications are just as commonplace with patients with CIEDs as the general population, therefore a full MR safety questionnaire should always be completed prior to entering the scan room.

### Consent

Informed written patient consent should be obtained and documented. However, local sites may decide that written consent is not necessary for lower risk scenarios as listed in [table 3](#). The process should include discussion of the potential risks of scanning based on the specific CIED characteristics, and the benefits of the scan. Information therefore must be available for the person taking consent (who may be the radiologist or cardiologist dependent on local protocols) regarding the patient, device, clinical necessity and feasible alternatives to MR scanning. Ideally, information about the risks should be provided to the patient in advance of the MRI via a written patient information leaflet. Suggested phrases for consent in more common scenarios are provided in [table 4](#). Other scenarios, not included in [table 4](#) and which depend on the CIED, patient and MRI characteristics can be viewed as a spectrum from lower to higher risk relative to these data ([table 3](#)).

### Device interrogation and programming

CIED programming changes need to be performed manually for MR Unlabelled generators, including disabling all advanced

features and all tachycardia therapies and shocks for defibrillators. This requires careful, individualised programming strategies that incorporate patient and CIED factors. For MR Unlabelled CIED programming, a monitor mode (OVO or ODO) should be used if there is an adequate underlying rhythm; or VOO, DOO if there is presence of bradycardia (<40 bpm).<sup>12 35</sup> It is important that the initial programmed CIED settings are recorded prior to programming for the MRI scan, in order to ensure appropriate settings are restored post-MRI and in the event of software reset.

### During scan

There should be at least one healthcare professional available within the MRI department for the duration of the scan who has basic life support training. Additionally, there should be personnel able to reprogramme the CIED if required in the department. In practice, this may be the same healthcare professional (commonly the cardiac physiologist or cardiologist). Personnel who are ACLS trained should be available in the hospital at the time of the scan.

Patients should be monitored verbally and with both continuous ECG and pulse oximetry monitoring for the duration they are in the MR Environment. When planning the MRI protocol, scans should be abbreviated and steps taken to reduce risk where feasible (eg, reduced SAR, choice of field strength), guided by the MR Safety Experts. In order to minimise the risk of needing to recall patients (and hence expose them to potential additional



clinical risk), there should be a mechanism in place to check completeness of the image acquisition in real-time prior to scan completion and removing the patient from the scanner.

#### Postscan

*MR Unlabelled* CIEDs should be interrogated and programmed back to pre-MRI initial settings. This is a more manual process than *MR Conditional* CIEDs and therefore requires careful checking of each parameter. If there is any significant change in a parameter, this should be communicated to the patient and ongoing follow-up should be arranged. Suggested lead parameters classified as significant are: a decrease in sensed P wave amplitude  $\geq 50\%$ ; a decrease in sensed R wave amplitude  $\geq 25\%$ ; an increase in capture threshold  $\geq 0.5$  V; an absolute change in pacing lead impedance  $\geq 50\Omega$ ; an absolute change in high-voltage lead impedance  $\geq 3\Omega$ ; a decrease in battery voltage  $\geq 0.04$  V. This is based on a small test-retest substudy of the MagnaSafe registry ( $n=30$ ). In this study, there were no P wave amplitude decreases  $\geq 1.0$  V, no R wave amplitude decreases  $\geq 2.0$  V and no pacing threshold increases of  $\geq 0.5$  mV—suggesting these sensitivity thresholds are real. Pacing lead impedance changes  $\geq 50\Omega$  were noted in 3.6% of leads; and shock impedance changes  $\geq 3\Omega$  were observed in 17.6% of defibrillator leads.<sup>35 40</sup>

#### OTHER CONSIDERATIONS

##### Emergency scanning of CIEDs

There are several medical conditions where emergency out of hours scanning may be requested for diagnosis and treatment planning (for example suspected spinal cord compression). However if such emergency scans are required in patients with active (functional, non-battery depleted) CIEDs, the same protocols must be followed as per elective scanning. There are no clinical circumstances where MRI without reprogramming of active CIEDs and adequate supervision can be recommended. Alternative imaging modalities should be explored, and for most conditions, treatment can be initiated empirically. There is recognition of need for prompt emergency MRI in some scenarios and the aim should be for provision of scans (at least regionally within a network) as soon as possible, although the standard pathways for device reprogramming and monitoring as detailed above should be followed.

##### Reporting suspected safety events

Previous work has highlighted an under-reporting safety events.<sup>50</sup> Any possible safety events related to the CIED should be reported to the Medicines and Healthcare products Regulatory Agency using the Yellow Card system with addition data collection for audit of local practice.<sup>77</sup>

##### Information for patients

Patients with *MR Conditional* devices have the right to expect access to MRI services locally where the clinical indication is reasonable. For patients with *MR Unlabelled* devices or *MR Conditional* devices where the conditions cannot be met, patients should encourage clinicians to collaborate with specialist centres to ensure access to MRI if there is clear potential benefit to mitigate the potential risk. Patients (and/or their carers) have a duty to make their referring clinician and the MRI department aware that they have a cardiac implantable electronic device to facilitate safe planning of the MRI scan. Because of the complexity of the service, this usually requires coordination of different medical teams and may introduce some delays, but this should

account for the clinical urgency of the scan. Where local services are not available (particularly for patients with *MR Unlabelled* devices), further information on specialist regional centres can be found at [www.mrimypacemaker.com](http://www.mrimypacemaker.com). Charitable and patient bodies including the Arrhythmia Alliance ([www.heartrhythmalliance.org](http://www.heartrhythmalliance.org)) and Cardiomyopathy UK ([www.cardiomyopathy.org](http://www.cardiomyopathy.org)) can provide support and guidance.

#### Areas for further research

Significant progress has been made in recent years to develop strict protocols for patients with CIEDs undergoing MRI. There are still significant logistical burdens associated with performing MRI for this patient group, and development of protocols and tools to alleviate these burdens are needed. There is a growing appreciation of a spectrum of risk associated with MRI for patients with *MR Unlabelled* CIEDs or *MR Conditional* CIEDs with unfulfilled MR conditions. A growing body of experience will help to inform clinical decision-making in individual scenarios. It is likely some scenarios will be regarded as having similar risk profiles to fully *MR Conditional* CIEDs, while other scenarios will have higher-risk profiles. Given that patients with *MR Unlabelled* CIEDs or *MR Conditional* CIEDs with unfulfilled MR conditions will require urgent diagnoses reliant on MRI for many decades, service provision should be developed to reduce the health inequality in MRI access. This will require appropriate design of infrastructure and health economic data to inform policy. Patients with CIEDs are typically not included in research trials which incorporate MRI, even if the CIED is *MR Conditional*. This compounds the health inequality that patients experience. MRI protocols should be developed for clinical trials and ethics submissions should reflect the changing practice for patients with CIEDs, particularly in the context of cancer or neurological disease.

#### CONCLUSION

These multisocietal Consensus Recommendations outline guidance for safe delivery of MRI to patients with CIEDs, and aim to improve provision to address current inequities of service delivery. The majority of patients with CIEDs should now be able to undergo MRI, although it is anticipated that patients with *MR Unlabelled* devices or *MR Conditional* devices where it is not possible to meet all the specified conditions a local risk-benefit decision is needed. Collaborative interdisciplinary working is required to facilitate safe workflows and these guidelines incorporate recommendations from all stakeholders, including patients, to drive widespread adoption and encourage service expansion.

#### Author affiliations

<sup>1</sup>Department of Cardiology, Barts Heart Centre, Barts Health NHS Trust, London, UK

<sup>2</sup>Institute of Health Informatics, University College London, London, UK

<sup>3</sup>Medical Physics, Guy's & St Thomas' NHS Foundation Trust, London, UK

<sup>4</sup>Representative for the British Institute of Radiology, London, UK

<sup>5</sup>Department of Medical Physics and Bioengineering, NHS Highland, Inverness, UK

<sup>6</sup>Representative for Institute of Physics and Engineering in Medicine, York, UK

<sup>7</sup>Society of Radiographers, London, UK

<sup>8</sup>Department of Radiology, Royal Bournemouth Hospital, Bournemouth, UK

<sup>9</sup>Representative for British Association of MR Radiographers, Sheffield, UK

<sup>10</sup>Medicines and Healthcare Products Regulatory Agency, London, UK

<sup>11</sup>Arrhythmia Alliance & Atrial Fibrillation Association, Stratford upon Avon, UK

<sup>12</sup>Department of Radiology, Imperial College London, London, UK

<sup>13</sup>Representative for Royal College of Radiologists, London, UK

<sup>14</sup>Cardiology, Royal Cornwall Hospitals NHS Trust, Truro, UK

<sup>15</sup>Representative for British Heart Rhythm Society, Chipping Norton, UK

<sup>16</sup>Radiology, Glasgow Royal Infirmary, Glasgow, UK

<sup>17</sup>Representative of the British Society of Cardiovascular Imaging and British Society of Cardiovascular CT, London, UK

<sup>18</sup>Department of Cardiovascular Imaging, Barts Heart Centre, Barts Health NHS Trust, London, UK

<sup>19</sup>Institute of Cardiovascular Science, University College London, London, UK

<sup>20</sup>Representative of British Cardiovascular Society, London, UK

**Twitter** Anish Bhuvu @mrmpacemaker and Trudie Lobban @#tcasl

**Contributors** AB drafted the work. CM chaired the Working Group. All authors were involved in the Working Group, reviewed the work critically, distributed and collated feedback from their respective professional societies. Each professional society reviewed the manuscript and approved the final version. All authors gave final approval of the version to be published.

**Funding** This study was funded by British Heart Foundation (HFHF\_016).

**Competing interests** AB and CM are recipients of a British Heart Foundation Grant (BHF HFHF\_016) to make MRI available for patients with pacemakers. All other authors have no conflict of interest.

**Patient consent for publication** Not applicable.

**Ethics approval** Not applicable.

**Provenance and peer review** Commissioned; externally peer reviewed.

**Supplemental material** This content has been supplied by the author(s). It has not been vetted by BMJ Publishing Group Limited (BMJ) and may not have been peer-reviewed. Any opinions or recommendations discussed are solely those of the author(s) and are not endorsed by BMJ. BMJ disclaims all liability and responsibility arising from any reliance placed on the content. Where the content includes any translated material, BMJ does not warrant the accuracy and reliability of the translations (including but not limited to local regulations, clinical guidelines, terminology, drug names and drug dosages), and is not responsible for any error and/or omissions arising from translation and adaptation or otherwise.

#### ORCID iDs

Anish Bhuvu <http://orcid.org/0000-0001-7532-7815>

Charlotte Manisty <http://orcid.org/0000-0003-0245-7090>

#### REFERENCES

- Indik JH, Gimbel JR, Abe H. HRS expert consensus statement on magnetic resonance imaging and radiation exposure in patients with cardiovascular implantable electronic devices. *Heart Rhythm* 2017;17:14:e99–151.
- Brignole M, Auricchio A, Baron-Esquivias G. 2013 ESC guidelines on cardiac pacing and cardiac resynchronization therapy: the task force on cardiac pacing and resynchronization therapy of the European Society of cardiology (ESC). developed in collaboration with the European heart rhythm association (EHRA). *Eur Heart J* 2013;34:2281–329.
- Cruyteninck Y, Dubourg B, Michelin P, et al. Pacemakers and MRI: a protocol in line with international guidelines and Approved by the SFICV (French Society of cardiovascular imaging). *Diagn Interv Imaging* 2017;98:203–15.
- Sommer T, Bauer W, Fischbach K, et al. Mr imaging in patients with cardiac pacemakers and implantable cardioverter defibrillators. *Fortschr Röntgenstr* 2017;189:204–17.
- Verma A, Ha ACT, Dennie C, et al. Canadian heart rhythm Society and Canadian association of radiologists consensus statement on magnetic resonance imaging with cardiac implantable electronic devices. *Can J Cardiol* 2014;30:1131–41.
- Vigen KK, Reeder SB, et al. Recommendations for imaging patients with cardiac implantable electronic devices (CIEDs). *J Magn Reson Imaging* 2020;1:1–7.
- Lowe MD, Plummer CJ, Manisty CH, et al. Safe use of MRI in people with cardiac implantable electronic devices. *Heart* 2015;101:1950–3.
- Stühlinger M, Burri H, Vernooy K. EHRA consensus on prevention and management of interference due to medical procedures in patients with cardiac implantable electronic devices. *Europace* 2022;00:1–26.
- Muthalaly RG, Nerlekar N, Ge Y, et al. Mri in patients with cardiac implantable electronic devices. *Radiology* 2018;289:281–92.
- Operational Information for Commissioning NE. Diagnostic imaging dataset. National health service. Available: <https://did.hscic.gov.uk>
- Murgatroyd F, Linker N, Cunningham D. National audit of cardiac rhythm management devices 2013–2021. Available: <http://www.bhrs.com/files/files/Audit Reports/2013%202014-CRM National Clinical Audit Report.pdf>
- Bradshaw PJ, Stobie P, Knuiman MW, et al. Trends in the incidence and prevalence of cardiac pacemaker insertions in an ageing population. *Open Heart* 2014;1:e000177.
- Sakakibara Y, Mitsui T. Concerns about sources of electromagnetic interference in patients with pacemakers. *Jpn Heart J* 1999;40:737–43.
- Celentano E, Caccavo V, Santamaria M. Access to magnetic resonance imaging of patients with magnetic resonance-conditional pacemaker and implantable cardioverter defibrillator systems: results from the really ProMRI study. *EP Eur* 2017;20:1001–9.
- Bhuvu AN, Moralee R, Moon JC, et al. Making MRI available for patients with cardiac implantable electronic devices: growing need and barriers to change. *Eur Radiol* 2020;30:1378–84.
- Nazarian S, Roguin A, Zviman MM, et al. Clinical utility and safety of a protocol for noncardiac and cardiac magnetic resonance imaging of patients with permanent pacemakers and implantable-cardioverter defibrillators at 1.5 Tesla. *Circulation* 2006;114:1277–84.
- Nazarian S, Reynolds MR, Ryan MP, et al. Utilization and likelihood of radiologic diagnostic imaging in patients with implantable cardiac defibrillators. *J Magn Reson Imaging* 2016;43:115–27.
- Sabzevari K, Oldman J, Herrey AS, et al. Provision of magnetic resonance imaging for patients with 'MR-conditional' cardiac implantable electronic devices: an unmet clinical need. *Europace* 2017;19:425–31.
- Clinical Imaging Board and British Cardiovascular Society. Joint statement: MRI for pacemakers and implantable cardioverter-defibrillators, 2018.
- Society BHR. Standards for implantation and follow-up of cardiac rhythm management devices in adults, 2018. Available: <http://www.bhrs.com/files/files/Guidelines/180122-sp-BHRS Standards - Implantation and Follow-Up of CRM Devices in Adults.pdf>
- Pieri C, Bhuvu A, Moralee R, et al. Access to MRI for patients with cardiac pacemakers and implantable cardioverter defibrillators. *Open Heart* 2021;8:e001598.
- NICE. NICE Guidance - Prostate cancer: diagnosis and management. *BJU Int* 2019;124:26 <http://doi.wiley.com/>
- Medicines and Healthcare products Regulatory Agency. Safety guidelines for magnetic resonance imaging equipment in clinical use, 2021. Available: <https://www.gov.uk/government/publications/safety-guidelines-for-magnetic-resonance-imaging-equipment-in-clinical-use>
- Shinbane JS, Colletti PM, Shellock FG. Magnetic resonance imaging in patients with cardiac pacemakers: era of "MR Conditional" designs. *J Cardiovasc Magn Reson* 2011;13:63.
- Wilkoﬀ BL, Bello D, Taborsky M, et al. Magnetic resonance imaging in patients with a pacemaker system designed for the magnetic resonance environment. *Heart Rhythm* 2011;8:65–73.
- Rod Gimbel J, Bello D, Schmitt M, et al. Randomized trial of pacemaker and lead system for safe scanning at 1.5 Tesla. *Heart Rhythm* 2013;10:685–91.
- Gold MR, Sommer T, Schwitzer J, et al. Full-Body MRI in patients with an implantable cardioverter-defibrillator. *J Am Coll Cardiol* 2015;65:2581–8.
- Williamson BD, Gohn DC, Ramza BM, et al. Real-World Evaluation of Magnetic Resonance Imaging in Patients With a Magnetic Resonance Imaging Conditional Pacemaker System: Results of 4-Year Prospective Follow-Up in 2,629 Patients. *JACC Clin Electrophysiol* 2017;3:1231–9.
- Kusumoto FM, Schoenfeld MH, Wilkoﬀ BL, et al. 2017 Hrs expert consensus statement on cardiovascular implantable electronic device lead management and extraction. *Heart Rhythm* 2017;14:e503–51.
- Wazni O, Epstein LM, Carrillo RG, et al. Lead extraction in the contemporary setting: the lexicon study: an observational retrospective study of consecutive laser lead extractions. *J Am Coll Cardiol* 2010;55:579–86.
- Jones SO, Eckart RE, Albert CM, et al. Large, single-center, single-operator experience with transvenous lead extraction: outcomes and changing indications. *Heart Rhythm* 2008;5:520–5.
- Poole JE, Gleva MJ, Mela T, et al. Complication rates associated with pacemaker or implantable cardioverter-defibrillator generator replacements and upgrade procedures: results from the replace registry. *Circulation* 2010;122:1553–61.
- The MRI my pacemaker campaign. Available: <https://www.mrimpacemaker.com/> [Accessed 04 Jul 2022].
- Bhuvu AN, Moralee R, Brunner T. Evidence to support magnetic resonance conditional labelling of all pacemaker and defibrillator leads in patients with cardiac implantable electronic devices. *Eur Heart J* 2021;ehab350.
- Russo RJ, Costa HS, Silva PD, et al. Assessing the risks associated with MRI in patients with a pacemaker or defibrillator. *N Engl J Med* 2017;376:755–64.
- Savoure A, Mechulan A, Burban M, et al. The KorA pacemaker is safe and effective for magnetic resonance imaging. *Clin Med Insights Cardiol* 2015;9:85–90.
- Schwitzer J, Gold MR, Al Fagih A, et al. Image quality of cardiac magnetic resonance imaging in patients with an implantable cardioverter defibrillator system designed for the magnetic resonance imaging environment. *Circ Cardiovasc Imaging* 2016;9:1–11.
- Schwitzer J, Kanal E, Schmitt M, et al. Impact of the Advia MRI pacing system on the diagnostic quality of cardiac MR images and contraction patterns of cardiac muscle during scans: Advia MRI randomized clinical multicenter study results. *Heart Rhythm* 2013;10:864–72.
- Bhuvu AN, Kellman P, Graham A, et al. Clinical impact of cardiovascular magnetic resonance with optimized myocardial scar detection in patients with cardiac implantable devices. *Int J Cardiol* 2019;279:72–8.
- Cohen JD, Costa HS, Russo RJ. Determining the risks of magnetic resonance imaging at 1.5 Tesla for patients with pacemakers and implantable cardioverter defibrillators. *Am J Cardiol* 2012;110:1631–6.
- Nowak B, Hemmer W, Israel CW. Safety of asynchronous ventricular pacemaker stimulation. statement of the Working group cardiac pacemakers of the German cardiac Society. *Clin Res Cardiol* 2006;95:57–60.

- 42 Elliott PM, Anastasakis A, Authors/Task Force members, *et al.* 2014 ESC guidelines on diagnosis and management of hypertrophic cardiomyopathy: the task force for the diagnosis and management of hypertrophic cardiomyopathy of the European Society of cardiology (ESC). *Eur Heart J* 2014;35:2733–79.
- 43 Bardy GH, Lee KL, Mark DB, *et al.* Amiodarone or an implantable cardioverter-defibrillator for congestive heart failure. *N Engl J Med* 2005;352:225–37.
- 44 Hilbert S, Jahnke C, Loebe S. Cardiovascular magnetic resonance imaging in patients with cardiac implantable electronic devices: a device-dependent imaging strategy for improved image quality. *Eur Heart J - Cardiovasc Imaging* 2017;1–11.
- 45 Roguin A, Zviman MM, Meininger GR, *et al.* Modern pacemaker and implantable Cardioverter/Defibrillator systems can be magnetic resonance imaging safe. *Circulation* 2004;110:475–82.
- 46 Nazarian S, Beinart R, Halperin HR. Magnetic resonance imaging and implantable devices. *Circ Arrhythm Electrophysiol* 2013;6:419–28.
- 47 Nazarian S, Hansford R, Rahsepar AA, *et al.* Safety of magnetic resonance imaging in patients with cardiac devices. *N Engl J Med* 2017;377:2555–64.
- 48 Gupta SK, Ya'qoub L, Wimmer AP, *et al.* Safety and clinical impact of MRI in patients with Non-MRI-conditional cardiac devices. *Radiol Cardiothorac Imaging* 2020;2:e200086.
- 49 Shah AD, Morris MA, Hirsh DS, *et al.* Magnetic resonance imaging safety in nonconditional pacemaker and defibrillator recipients: a meta-analysis and systematic review. *Heart Rhythm* 2018;15:1001–8.
- 50 Irmich W, Irmich B, Bartsch C, *et al.* Do we need pacemakers resistant to magnetic resonance imaging? *Europace* 2005;7:353–65.
- 51 Munawar DA, Chan JEZ, Emami M, *et al.* Magnetic resonance imaging in non-conditional pacemakers and implantable cardioverter-defibrillators: a systematic review and meta-analysis. *Europace* 2020;22:288–98.
- 52 Wilkoff BL, Albert T, Lazebnik M, *et al.* Safe magnetic resonance imaging scanning of patients with cardiac rhythm devices: a role for computer modeling. *Heart Rhythm* 2013;10:1815–21 <https://linkinghub.elsevier.com/retrieve/pii/S1547527113011466> doi:10.1016/j.hrthm.2013.10.009
- 53 Gold MR, Kanal E, Schwittler J, *et al.* Preclinical evaluation of implantable cardioverter-defibrillator developed for magnetic resonance imaging use. *Heart Rhythm* 2015;12:631–8.
- 54 Higgins JV, Sheldon SH, Watson RE, *et al.* "Power-on resets" in cardiac implantable electronic devices during magnetic resonance imaging. *Heart Rhythm* 2015;12:540–4.
- 55 Al-Dayeh L, Rahman M, Venook R. Practical aspects of Mr imaging safety test methods for Mr conditional active implantable medical devices. *Magn Reson Imaging Clin N Am* 2020;28:559–71.
- 56 Achenbach S, Moshage W, Diem B, *et al.* Effects of magnetic resonance imaging on cardiac pacemakers and electrodes. *Am Heart J* 1997;134:467–73.
- 57 Wang Y, Guo R, Hu W, *et al.* Magnetic resonance conditionality of abandoned leads from active implantable medical devices at 1.5 T. *Magn Reson Med* 2022;87:394–408.
- 58 Luechinger R, Zeijlemaker VA, Pedersen EM, *et al.* In vivo heating of pacemaker leads during magnetic resonance imaging. *Eur Heart J* 2005;26:376–83.
- 59 Langman DA, Goldberg IB, Finn JP, *et al.* Pacemaker lead tip heating in abandoned and pacemaker-attached leads at 1.5 Tesla MRI. *J Magn Reson Imaging* 2011;33:426–31.
- 60 Mattei E, Gentili G, Censi F, *et al.* Impact of capped and uncapped abandoned leads on the heating of an MR-conditional pacemaker implant. *Magn Reson Med* 2015;73:390–400.
- 61 Golestanirad L, Rahsepar AA, Kirsch JE, *et al.* Changes in the specific absorption rate (SAR) of radiofrequency energy in patients with retained cardiac leads during MRI at 1.5T and 3T. *Magn Reson Med* 2019;81:653–69.
- 62 Yao A, Goren T, Samaras T, *et al.* Radiofrequency-induced heating of broken and abandoned implant leads during magnetic resonance examinations. *Magn Reson Med* 2021;86:2156–2164.
- 63 Padmanabhan D, Kella DK, Mehta R, *et al.* Safety of magnetic resonance imaging in patients with legacy pacemakers and defibrillators and abandoned leads. *Heart Rhythm* 2018;15:228–33.
- 64 Higgins JV, Gard JJ, Sheldon SH, *et al.* Safety and outcomes of magnetic resonance imaging in patients with abandoned pacemaker and defibrillator leads. *Pacing Clin Electrophysiol* 2014;37:1284–90.
- 65 Schaller RD, Brunner T, Riley MP, *et al.* Magnetic resonance imaging in patients with cardiac implantable electronic devices with abandoned leads. *JAMA Cardiol* 2021;6:549–8.
- 66 Hartnell GG, Spence L, Hughes LA, *et al.* Safety of Mr imaging in patients who have retained metallic materials after cardiac surgery. *AJR Am J Roentgenol* 1997;168:1157–9.
- 67 Pulver AF, Puchalski MD, Bradley DJ, *et al.* Safety and imaging quality of MRI in pediatric and adult congenital heart disease patients with pacemakers. *Pacing Clin Electrophysiol* 2009;32:450–6.
- 68 Kanal E. Safety of Mr imaging in patients with retained epicardial PACER wires. *AJR Am J Roentgenol* 1998;170:213–4.
- 69 Gatterer C, Stelzmüller M-E, Kammerlander A, *et al.* Safety and image quality of cardiovascular magnetic resonance imaging in patients with retained epicardial pacing wires after heart transplantation. *J Cardiovasc Magn Reson* 2021;23:24.
- 70 Vuorinen A-M, Pakarinen S, Jaakkola I, *et al.* Clinical experience of magnetic resonance imaging in patients with cardiac pacing devices: unrestricted patient population. *Acta Radiol* 2019;60:1414–21.
- 71 Vuorinen A-M, Lehtonen J, Pakarinen S, *et al.* Cardiac magnetic resonance imaging-based screening for cardiac sarcoidosis in patients with atrioventricular block requiring temporary pacing. *J Am Heart Assoc* 2022;11:e024257.
- 72 Friedman HL, Acker N, Dalzell C, *et al.* Magnetic resonance imaging in patients with recently implanted pacemakers. *Pacing Clin Electrophysiol* 2013;36:1090–5.
- 73 Sommer T, Naehle CP, Yang A, *et al.* Strategy for safe performance of extrathoracic magnetic resonance imaging at 1.5 tesla in the presence of cardiac pacemakers in non-pacemaker-dependent patients: a prospective study with 115 examinations. *Circulation* 2006;114:1285–92.
- 74 Okamura H, Padmanabhan D, Watson RE. Magnetic resonance imaging in non-dependent pacemaker patients with pacemakers and defibrillators with a nearly depleted battery. *Pacing Clin Electrophysiol* 2017;1–29.
- 75 Rahsepar AA, Zimmerman SL, Hansford R, *et al.* The relationship between MRI radiofrequency energy and function of nonconditional implanted cardiac devices: a prospective evaluation. *Radiology* 2020;295:307–13.
- 76 Mollerus M, Albin G, Lipinski M, *et al.* Magnetic resonance imaging of pacemakers and implantable cardioverter-defibrillators without specific absorption rate restrictions. *Europace* 2010;12:947–51.
- 77 Yellow Card Scheme. Medicines & Healthcare Products Regulatory Agency. Available: <https://yellowcard.mhra.gov.uk/information> [Accessed 04 Jul 2022].
- 78 Nyotowidjojo IS, Skinner K, Shah AS, *et al.* Thoracic versus nonthoracic MR imaging for patients with an Mr nonconditional cardiac implantable electronic device. *Pacing Clin Electrophysiol* 2018;41:589–96.
- 79 Awad K, Griffin J, Crawford TC, *et al.* Clinical safety of the Iforia implantable cardioverter-defibrillator system in patients subjected to thoracic spine and cardiac 1.5-T magnetic resonance imaging scanning conditions. *Heart Rhythm* 2015;12:2155–61.
- 80 Naehle CP, Meyer C, Thomas D, *et al.* Safety of brain 3-T MR imaging with transmit-receive head coil in patients with cardiac pacemakers: pilot prospective study with 51 examinations. *Radiology* 2008;249:991–1001.
- 81 Gimbel JR. Magnetic resonance imaging of implantable cardiac rhythm devices at 3.0 Tesla. *Pacing Clin Electrophysiol* 2008;31:795–801.